

Omeprazole on Lux 5µ Cellulose-2 in RP Elution Mode (0.1% FA)

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4457-E0

Elution Type: Isocratic

Eluent A: 0.1% FA in H₂O

Eluent B: 0.1% FA in CH₃CN

Gradient

Step No.	Time (min)	Pct A	Pct B
1	20	70	30

Profile:

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20713



Products used in this application:

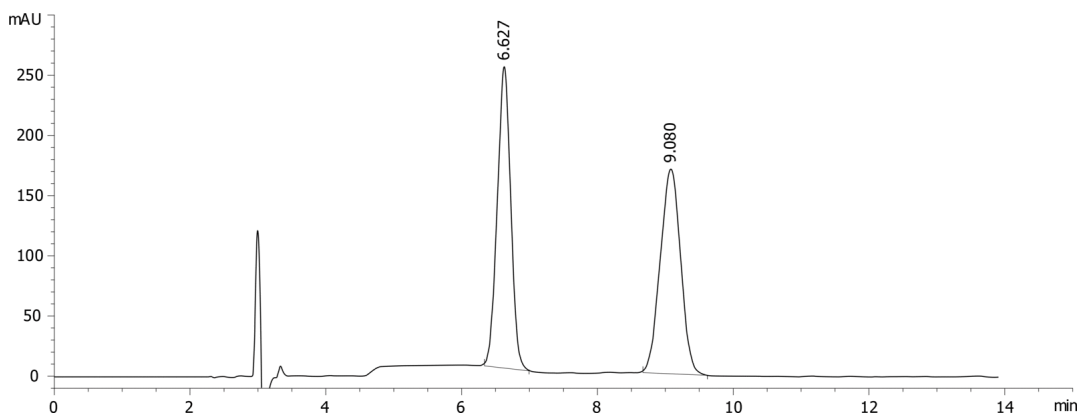

ANALYTES:

1 Omeprazole

2 Omeprazole

