

Tenatoprazole on Lux 5µm Cellulose-2 in RP Elution Mode (20mM NH₄HCO₃)

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4457-E0
Elution Type: Isocratic
Eluent A: 20mM NH₄HCO₃
Eluent B: 0.1%FA in CH₃CN
Gradient


Step No.	Time (min)	Pct A	Pct B
1	20	40	60

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20721



Products used in this application:

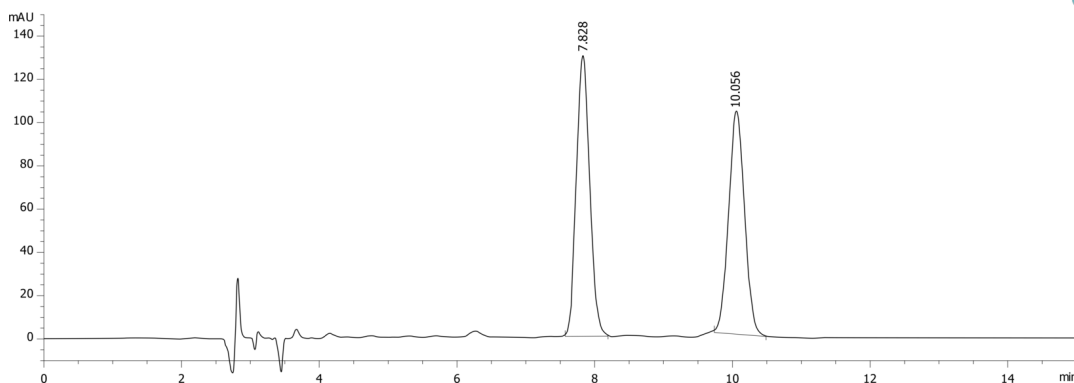

ANALYTES:

- 1 Tenatoprazole
- 2 Tenatoprazole

