

Tenatoprazole- lux 5µ Cellulose-2 in RP Elution Mode (0.1% FA)

Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4457-E0

Elution Type: Isocratic

Eluent A: 0.1% FA in H₂O

Eluent B: 0.1%FA in CH₃CN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	20	40	60

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

20723



Products used in this application:

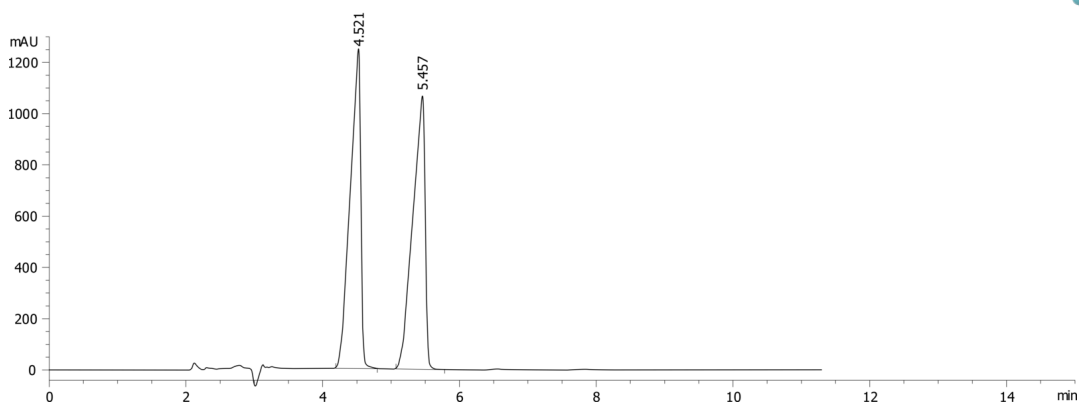

ANALYTES:

- 1 Tenatoprazole
- 2 Tenatoprazole

