

Prodigy Test Mix on Kinetex 1.7u C18 on Acquity System

Column: Kinetex® 1.7 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea

Dimensions: 50 x 2.1 mm ID

Order No: 00B-4475-AN

Elution Type: Isocratic

Eluent A: Water


Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	50	50

Flow Rate: 0.5 mL/min

Col. Temp.: 30 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ (ambient)

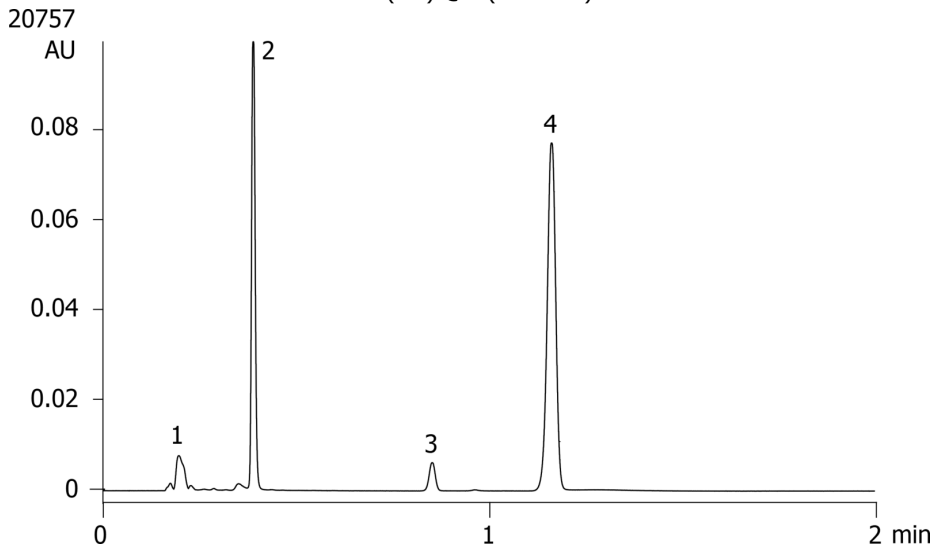

ANALYTES:

1 Uracil

Retention Time: 0.203 min

2 Acetophenone

Retention Time: 0.395 min

3 Toluene

Retention Time: 0.857 min

4 Naphthalene

Retention Time: 1.165 min



Products used in this application:

