

Acetylsalicylic acid using Kinetex 5u C18 250x4.6mm

Column: Kinetex® 5µm C18 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4601-E0

Elution Type: Isocratic

Eluent A: 600/400/2 Water/Acetonitrile/H3PO4

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 237 nm (ambient)

20770



Products used in this application:

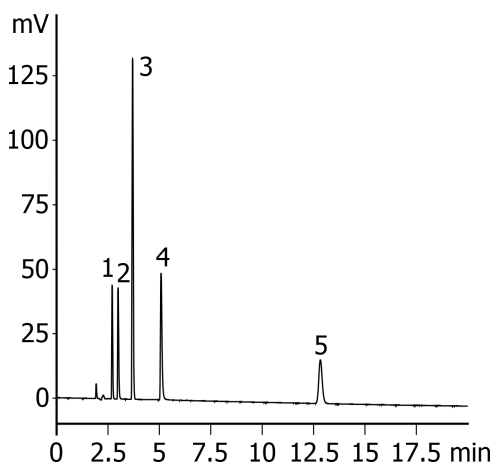

ANALYTES:

- 1 Impurity A
- 2 Impurity B
- 3 Impurity C
- 4 Acetylsalicylic acid
- 5 Impurity D

