

Preparative loading on Kinetex 5u C18 in Axia 50x21.2mm

Column: Kinetex® 5µm C18 100 Å, AXIA Packed LC Column 50 x 21.2 mm, Ea

Dimensions: 50 x 21.2 mm ID

Order No: 00B-4601-P0-AX

Elution Type: Gradient

Eluent A: Water with 0.5% Formic Acid

Eluent B: Acetonitrile with 0.5% Formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	8	50	50
	3	11	50	50

Flow Rate: 30 mL/min


Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

20774



Products used in this application:

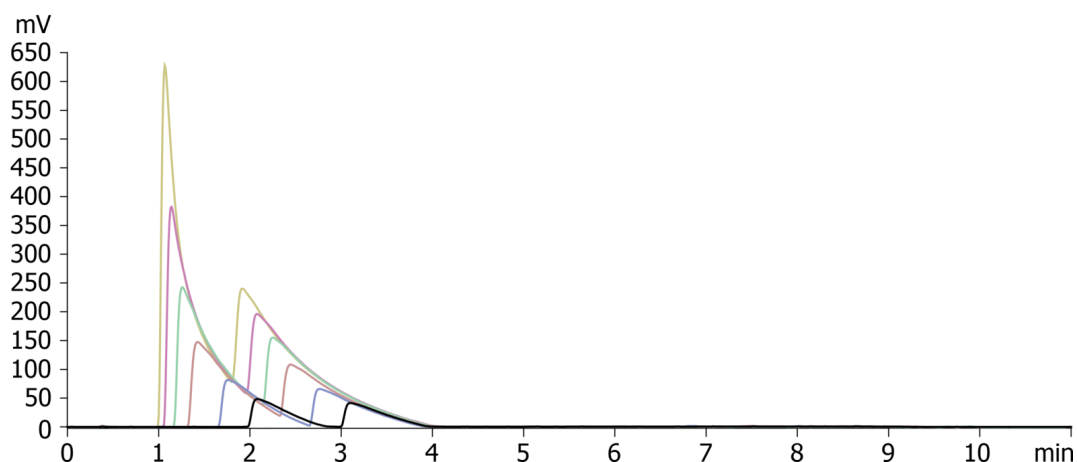

ANALYTES:

- 1 Doxepin
- 2 Amitriptyline

