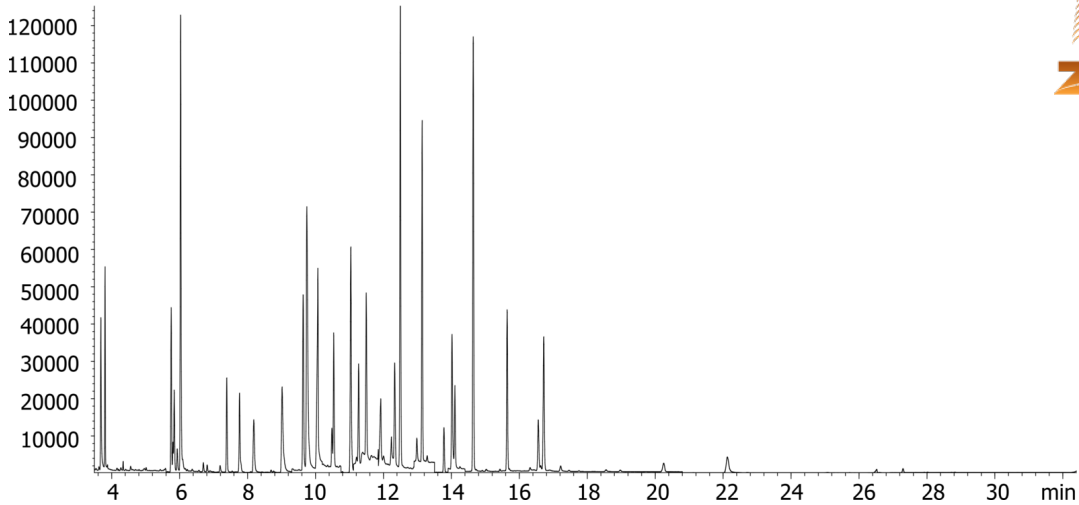


Quechers [AOAC-GCMS] - Lettuce extract, ZB-5MSi 30x0.25x0.25mm**Column:** Zebtron™ ZB-5MSi w/GUARDIAN™, 5 m, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl 95% Dimethylpolysiloxane**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G018-11-GGA**Oven Profile:** Initial temp: 100 C
Initial Hold Time: 0 min**Carrier Gas:** Constant Flow Helium, 1 mL/min**Injection:** Splitless :1 2 µL @ 250°C**Detection:** Mass Spectrometer (MS) (350°C)

120826

**Products used in this application:**

Quechers [AOAC-GCMS] - Lettuce extract, ZB-5MSi 30x0.25x0.25mm

ANALYTES:

- 1 Dichlorvos
- 2 o-Phenylphenol
- 3 Trifluralin
- 4 alpha-HCH-d6
- 5 Atrazine
- 6 Chlorothalonil
- 7 Chlorpyrifos-methyl
- 8 Carbaryl
- 9 Parathion-d10
- 10 Chlorpyrifos
- 11 Cyprodinil
- 12 Tolyfluanid
- 13 Procymidone
- 14 o,p-DDD
- 15 Kresoxim-methyl
- 16 Ethion
- 17 Endosulfan sulfate
- 18 Tebuconazole
- 19 Triphenyl phosphate
- 20 Bifenthrin
- 21 L-Cyhalothrin
- 22 Permethrins
- 23 Azoxystrobin

