

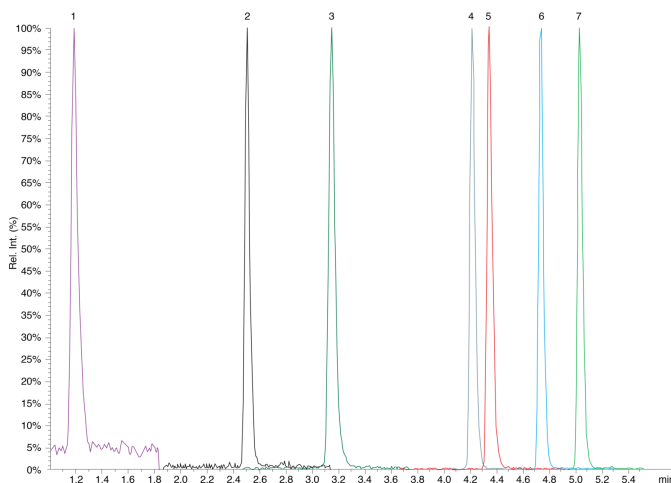
Quantification of Triazine Pesticides in Water by LC/MS/MS**Column:** Kinetex® 2.6 µm C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4462-AN**Elution Type:** Gradient**Eluent A:** Water with 0.1 % Formic Acid**Eluent B:** 100% MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	0.25	60	40
	3	2	60	40
	4	3	25	75
	5	4	25	75
	6	4.1	95	5

Flow Rate: 0.3 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (600 °C)**Detector Info:** <a target="_blank"

20883

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

**Products used in this application:****ANALYTES:****1** Atrazine-desethyl-desisopropyl

Retention Time: 1.19 min

2 Atrazine-desisopropyl

Retention Time: 2.5 min

3 Atrazine-desethyl

Retention Time: 3.1 min

4 Cyanazine

Retention Time: 4.22 min

5 Simazine

Retention Time: 4.34 min

6 Atrazine

Retention Time: 4.74 min

7 Propazine

Retention Time: 5.03 min

