


Chlorinated Pesticides (508 list) on ZB-SemiVolatiles**Column:** Zebron™ ZB-SemiVolatiles, GC Cap. Column 30 m x 0.25 mm x 0.25 µm, Ea**Phase:** 5% Phenyl-Arylene 95% Dimethylpolysiloxane (Enviro-Inert)**Dimensions:** 30 meters x 0.25 mm x 0.25 µm**Order No:** 7HG-G027-11**Oven Profile:** 80 °C for 0.75 min to 190 °C @ 35 °C/min to 240 °C @ 5 °C/min to 300 °C @ 20 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 1.2 mL/min**Injection:** Pulsed :1 1 µL @ 260°C**Detection:** Mass Selective (MSD) (320°C)**Analyst Note:** Pulsed Splitless Injection @ 30 psi for 0.55 min
Liner: Single Taper with Wool @ Bottom (AG0-8499)
Gold Seal: Easy Seal (AG0-8620)
Septum: PhenoRed (AG0-4696)**Products used in this application:**

Chlorinated Pesticides (508 list) on ZB-SemiVolatiles**ANALYTES:**

- 1** Hexachlorocyclopentadiene
- 2** Etridiazole
- 3** Chloroneb
- 4** Propachlor
- 5** Trifluralin
- 6** alpha-BHC
- 7** Hexachlorobenzene
- 8** Simazine
- 9** Atrazine
- 10** gamma-BHC
- 11** Pentachloronitrobenzene
- 12** b-BHC
- 13** Chlorothalonil
- 14** delta-BHC
- 15** Metribuzin
- 16** Alachlor
- 17** Heptachlor
- 18** Metolachlor
- 19** Cyanazine
- 20** Dacthal
- 21** Aldrin
- 22** 4,4-Dibromophenol
- 23** Heptachlor epoxide
- 24** gamma-Chlordane
- 25** alpha-Chlordane
- 26** Endosulfan II
- 27** DDE
- 28** Dieldrin
- 29** Endrin
- 30** Chlorobenzilate
- 31** Endosulfan II
- 32** DDD
- 33** Endrin aldehyde
- 34** Endosulfan sulfate
- 35** DDT
- 36** Methoxychlor
- 37** cis-Permethrin
- 38** trans-Permethrin

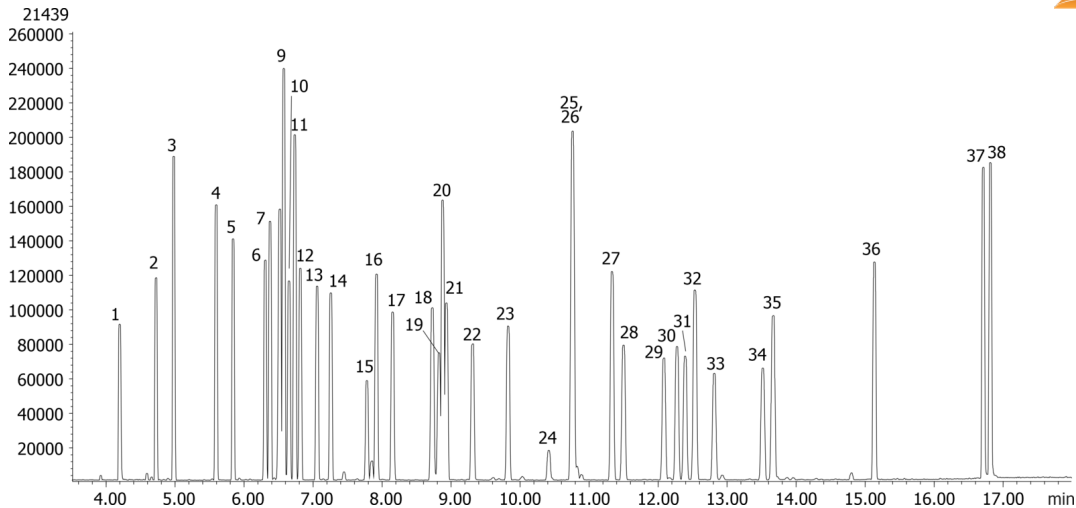