

Kinetex 5u XB-C18 Axia Efficiency Comparison

Column: Kinetex® 5µm XB-C18 100 Å, AXIA Packed LC Column 150 x 21.2 mm, Ea

Dimensions: 150 x 21.2 mm ID

Order No: 00F-4605-P0-AX

Elution Type: Isocratic

Eluent A: Acetonitrile

Eluent B: Water

Gradient **Step No.** **Time (min)** **Pct A** **Pct B**

Profile: **1** 0 50 50


Flow Rate: 25 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Products used in this application:

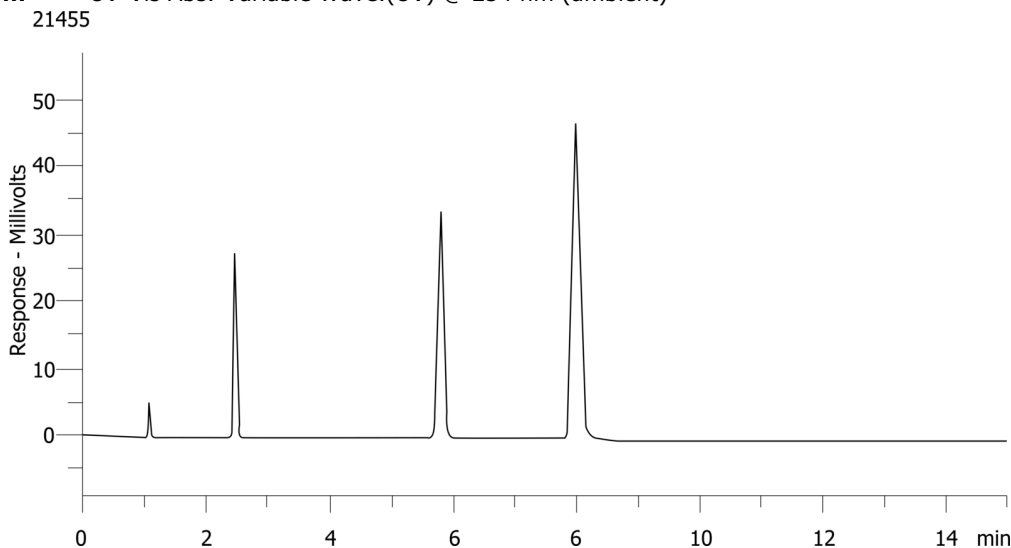

ANALYTES:

- 1 Uracil
- 2 Acetophenone
- 3 Toluene
- 4 Naphthalene

