


Yarra 3µm SEC-3000 Different Salt Concentrations: Lysozyme

Column: Yarra™ 3 µm SEC-3000, LC Column 300 x 7.8 mm, Ea
Dimensions: 300 x 7.8 mm ID
Order No: 00H-4513-K0
Elution Type: Isocratic
Eluent A: 50 mM Sodium Phosphate pH 6.8 / 0.3 M Sodium Chloride
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

21463

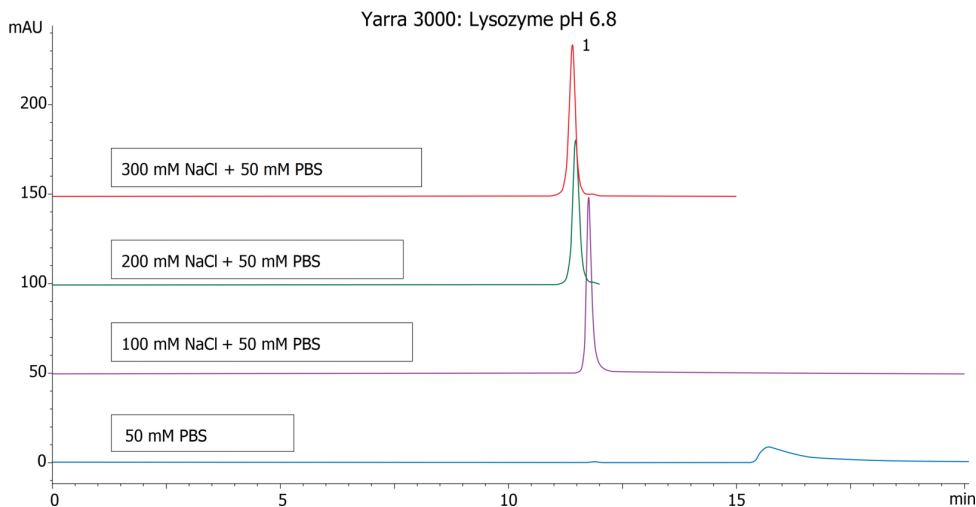

ANALYTES:

1 lysozyme

Retention Time: 11.5 min



Products used in this application:

