

## PVP on 2 X Phenogel 5 µm Linear(2) 300 X 7.8

**Column:** Phenogel™ 5 µm Linear(2), LC Column 300 x 7.8 mm, Ea

**Dimensions:** 300 x 7.8 mm ID

**Order No:** 00H-3259-K0

**Elution Type:** Isocratic

**Eluent A:** 10 mM LiBr in DMF

| Gradient | Step No. | Time (min) | Pct A |
|----------|----------|------------|-------|
| Profile: | 1        | 0          | 100   |


**Flow Rate:** 2 mL/min

**Col. Temp.:** 40 °C

**Detection:** Refractive Index (RI) @ (40 °C)  
21505



Products used in this application:

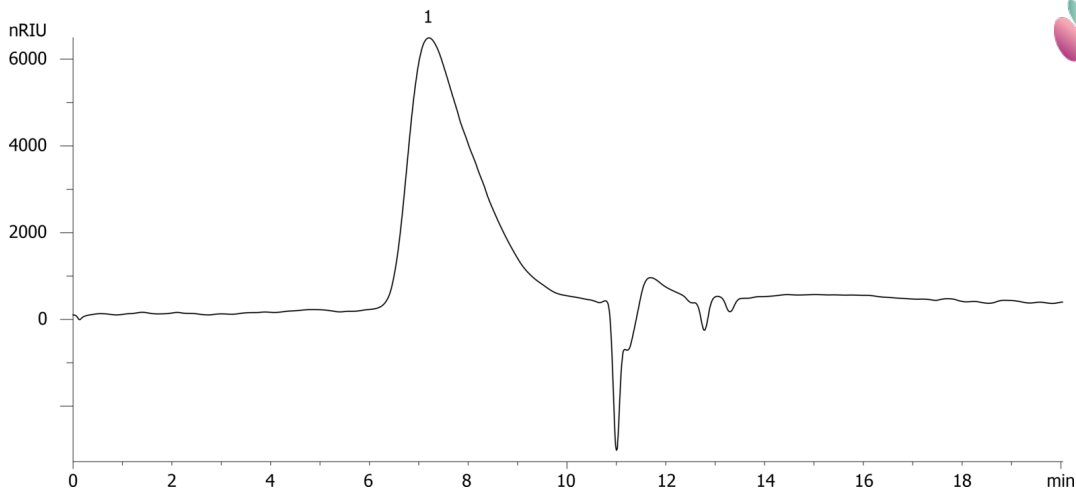

### ANALYTES:

1 Polyvinylpyrrolidone

Retention Time: 7.2 min

