

FMOC-Val-OH on Lux 5µm Cellulose-3 in NP

Column: Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4493-E0

Elution Type: Isocratic

Eluent A: Hexane / EtOH (80:20) Formic acid (0.1%)


Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

21535

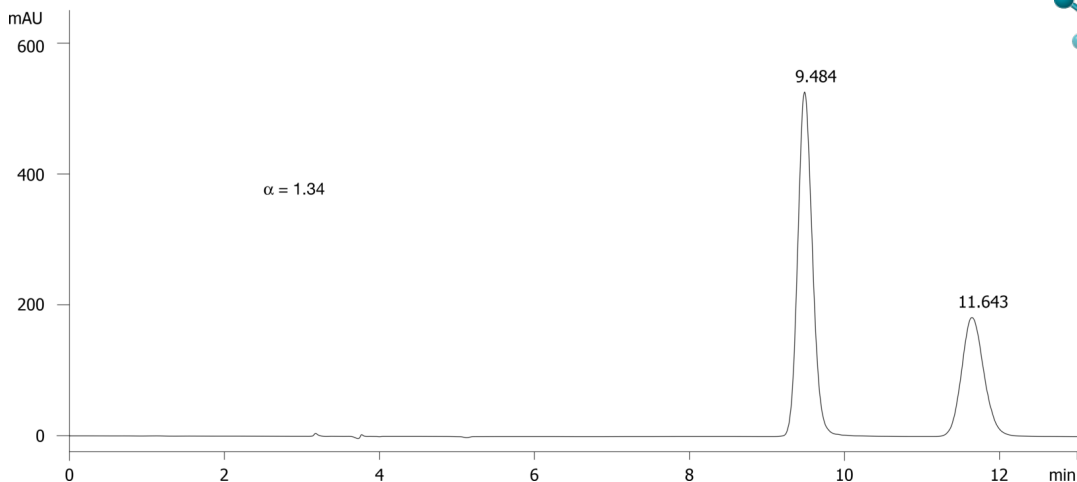

Products used in this application:



ANALYTES:

1 FMOC-Val-OH

Retention Time: 9.484 min

2 FMOC-Val-OH

Retention Time: 11.643 min

