

FMOC-Met-OH on Lux 5µm Cellulose-3 in NP


Column: Lux® 5 µm Cellulose-3, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4493-E0
Elution Type: Isocratic
Eluent A: Hexane / EtOH (80:20) Formic acid (0.1%)
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)



Products used in this application:

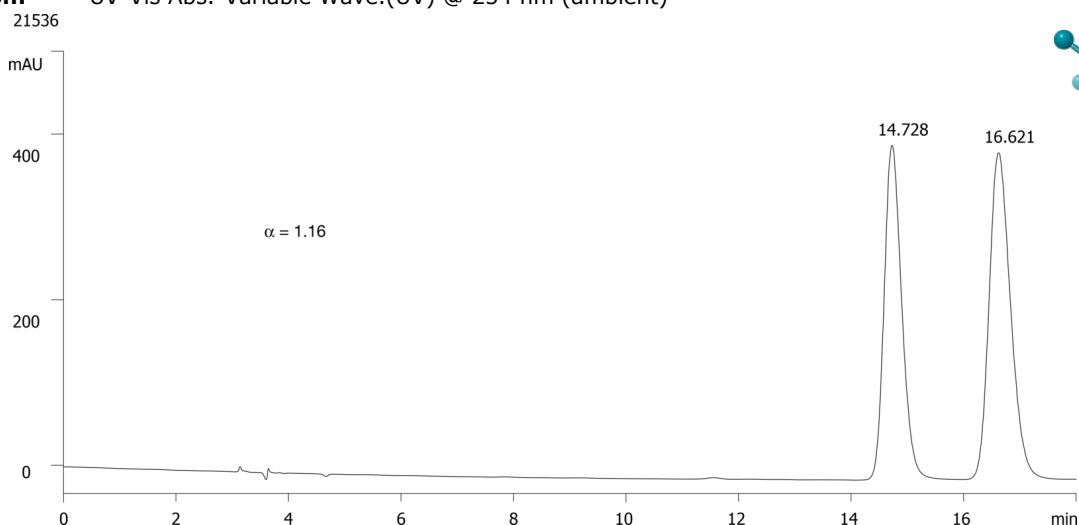

ANALYTES:

- 1 FMOC-Met-OH
Retention Time: 14.728 min
- 2 FMOC-Met-OH
Retention Time: 16.621 min

