

FMOC-Glu(OtBu)-OH on Lux 5 μ m Cellulose-1 in RP

Column: Lux[®] 5 μ m Cellulose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4459-E0

Elution Type: Isocratic

Eluent A: ACN / 0.1% TFA (60:40) with TFA (0.1%) in buffer

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

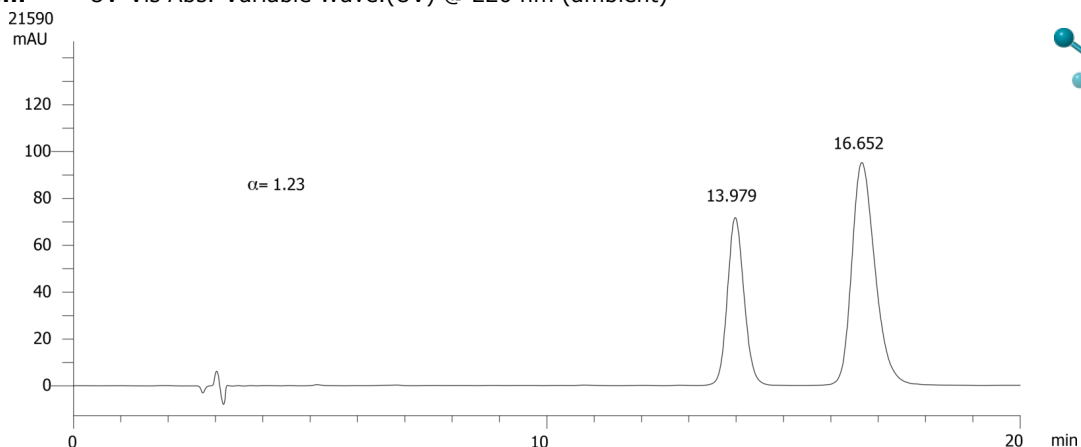

ANALYTES:

- 1 FMOC-Glu(OtBu)-OH
Retention Time: 13.979 min
- 2 FMOC-Glu(OtBu)-OH
Retention Time: 16.652 min

