

Epoxiconazole on Lux 5µm Cellulose-4 in NP

Column: Lux® 5 µm Cellulose-4, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4491-E0

Elution Type: Isocratic

Eluent A: Hexane / IPA (50:50) DEA (0.1%)

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

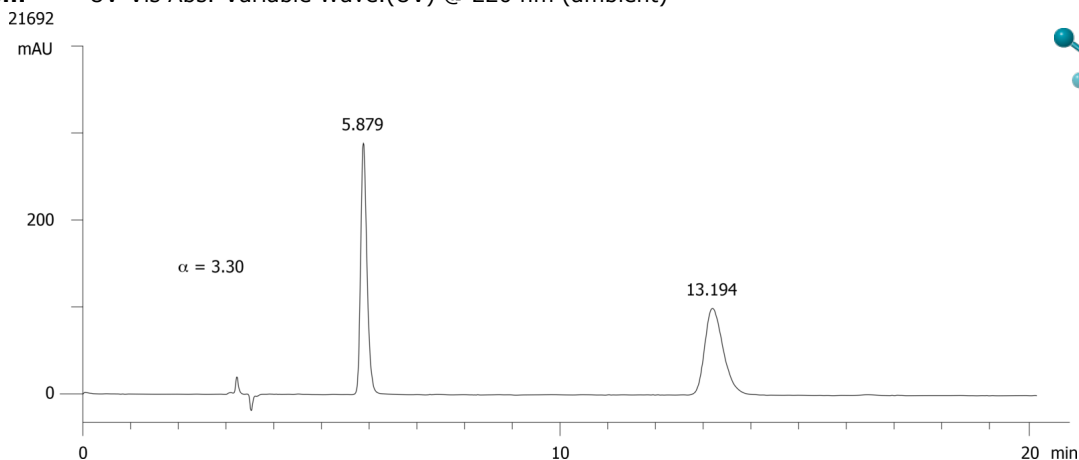

ANALYTES:

1 Epoxiconazole

Retention Time: 5.879 min

2 Epoxiconazole

Retention Time: 13.194 min

