

Bromuconazole on Lux 5µm Cellulose-4 in RP


Column: Lux® 5 µm Cellulose-4, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4491-E0
Elution Type: Isocratic
Eluent A: ACN / 20mM NH₄HCO₃+0.1%DEA (60:40) with DEA (0.1%) in buffer
Gradient

Step No.	Time (min)	Pct A
1	0	100

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)
 App ID 21752



Products used in this application:

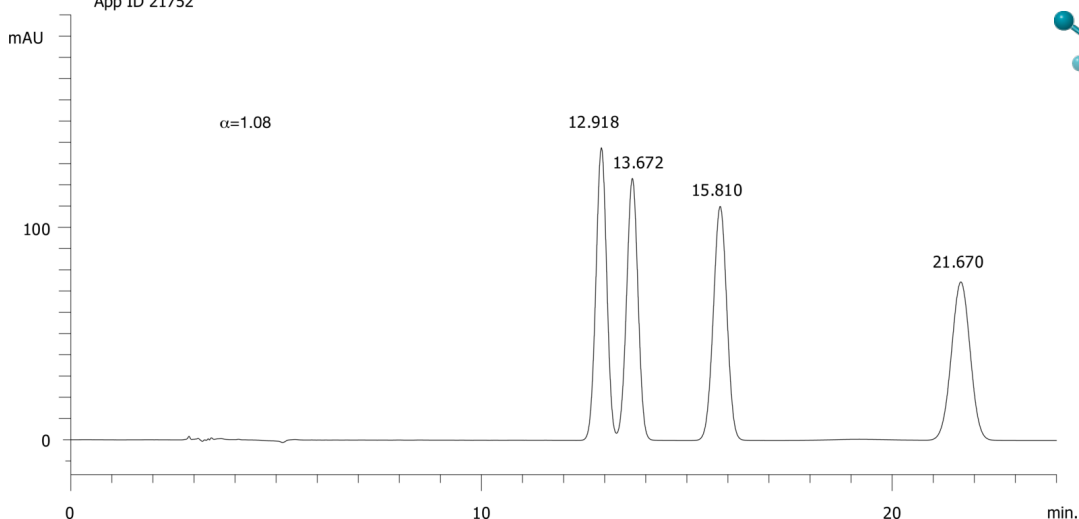

ANALYTES:

- 1** Bromuconazole
Retention Time: 12.918 min
- 2** Bromuconazole
Retention Time: 13.672 min
- 3** Bromuconazole
Retention Time: 15.81 min
- 4** Bromuconazole
Retention Time: 21.67 min

