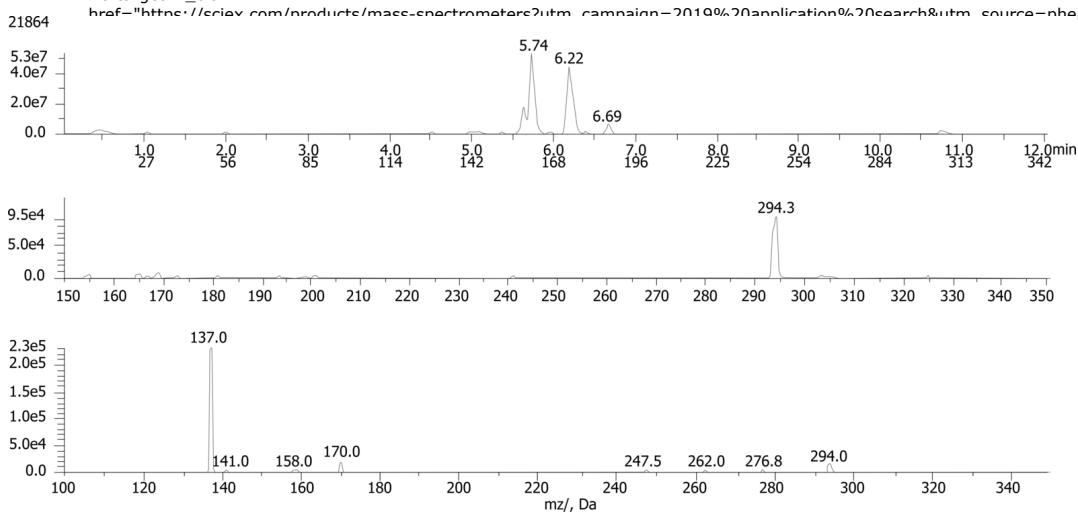


Precursor Scan of Nordihydrocapsaicin LC/MS/MS on Kinetex 5u C18 100x2.1**Column:** Kinetex® 5µm C18 100 Å, LC Column 100 x 2.1 mm, Ea**Dimensions:** 100 x 2.1 mm ID**Order No:** 00D-4601-AN**Elution Type:** Gradient**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 0.1% Formic Acid in Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	70	30
	2	7	10	90
	3	9	10	90
	4	9.1	70	30
	5	12	70	30

Flow Rate: 0.6 mL/min**Col. Temp.:** 25 °C**Detection:** Mass Spectrometer (MS) @ amu (600 °C)**Detector Info:** SCIEX<**Products used in this application:****ANALYTES:****1** Nordihydrocapsaicin