

Fmoc-Asn(Trt)-OH on Lux Cellulose-2 250x4.6mm in RP


Column: Lux® 5 µm Cellulose-2, LC Column 250 x 4.6 mm, Ea
Dimensions: 250 x 4.6 mm ID
Order No: 00G-4457-E0
Elution Type: Isocratic
Eluent A: Water + 0.1%TFA
Eluent B: Acetonitrile + 0.1% TFA
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	45	55

Profile:
Flow Rate: 1 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Variable Wave.(UV) @ (ambient) 21873



Products used in this application:

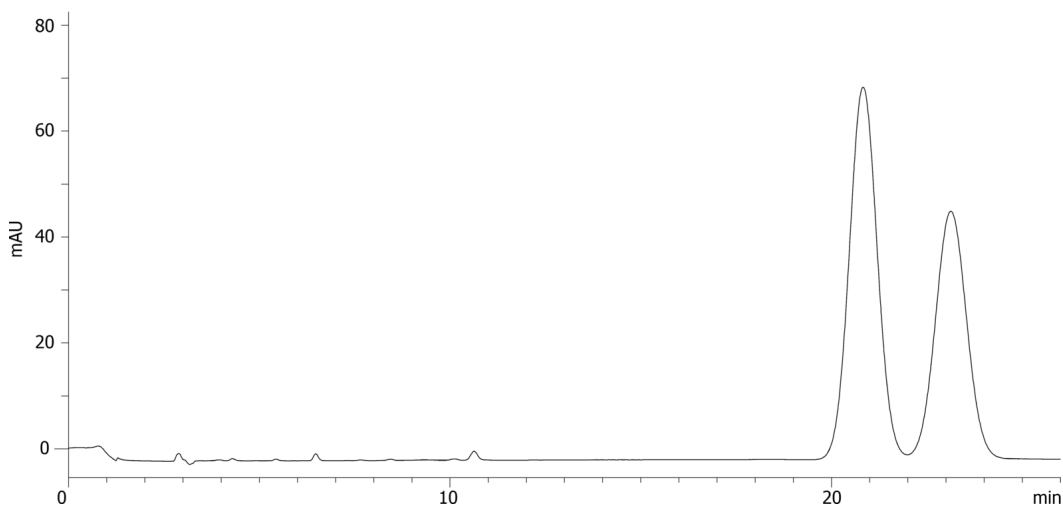

ANALYTES:

- 1** FMOC-Asn(Trt)-OH
Retention Time: 20.825 min
- 2** FMOC-Asn(Trt)-OH
Retention Time: 23.124 min

