

Paroxetine EP method on Kinetex 5µm C8 250x4.6 mm

Column: Kinetex® 5µm C8 100 Å, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4608-E0

Elution Type: Gradient

Eluent A: TFA/THF/Water (5:100:900)

Eluent B: TFA/THF/ACN (5:100:900)

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	30	80	20
	3	50	20	80
	4	55	20	80
	5	60	80	20
	6	65	80	20

Flow Rate: 1 mL/min


Col. Temp.: 40 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 295 nm (ambient)

21927



Products used in this application:

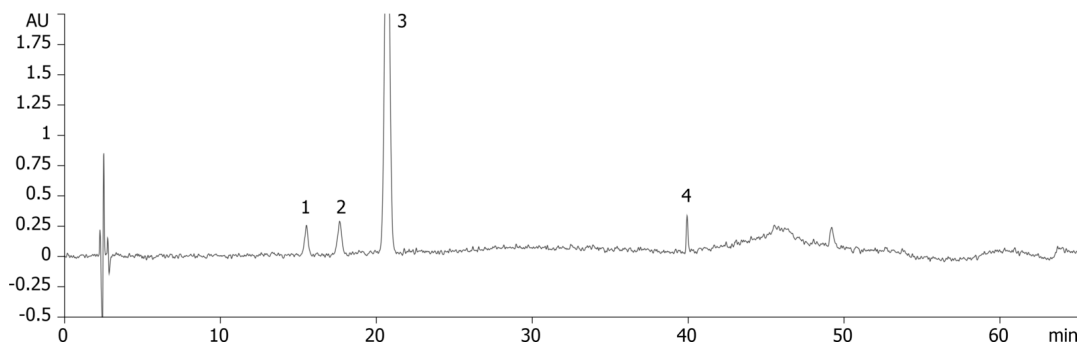

ANALYTES:

- 1 Impurity A
- 2 Impurity E
- 3 Paroxetine
- 4 Impurity C

