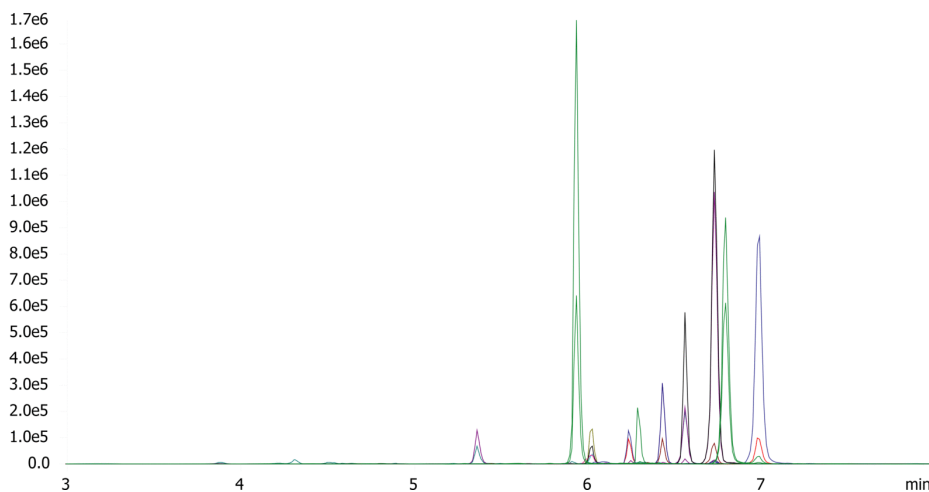


Mycotoxin analysis using Kinetex 2.6u Biphenyl 50x2.1mm**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4622-AN**Elution Type:** Gradient**Eluent A:** 5mM ammonium acetate with 0.1% acetic acid**Eluent B:** Methanol with 5mM ammonium acetate with 0.1% acetic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2
	2	2	98	2
	3	5	20	80
	4	5.2	2	98
	5	8	2	98
	6	8.01	98	2
	7	11	98	2

Flow Rate: 0.45 mL/min**Col. Temp.:** ambient**Detection:** Mass Spectrometer (MS) @ amu (ambient)**Detector Info:** SCIEX
21951**Products used in this application:**

Mycotoxin analysis using Kinetex 2.6u Biphenyl 50x2.1mm

ANALYTES:

- 1** 15-Acetyldeoxynivalenol
- 2** DAS
- 3** FB1
- 4** HT2 Toxin
- 5** FB2
- 6** T2 Toxin
- 7** Aflatoxin M1
- 8** Aflatoxin G2
- 9** Ochratoxin A
- 10** Aflatoxin G1
- 11** Aflatoxin B2
- 12** Aflatoxin B1

