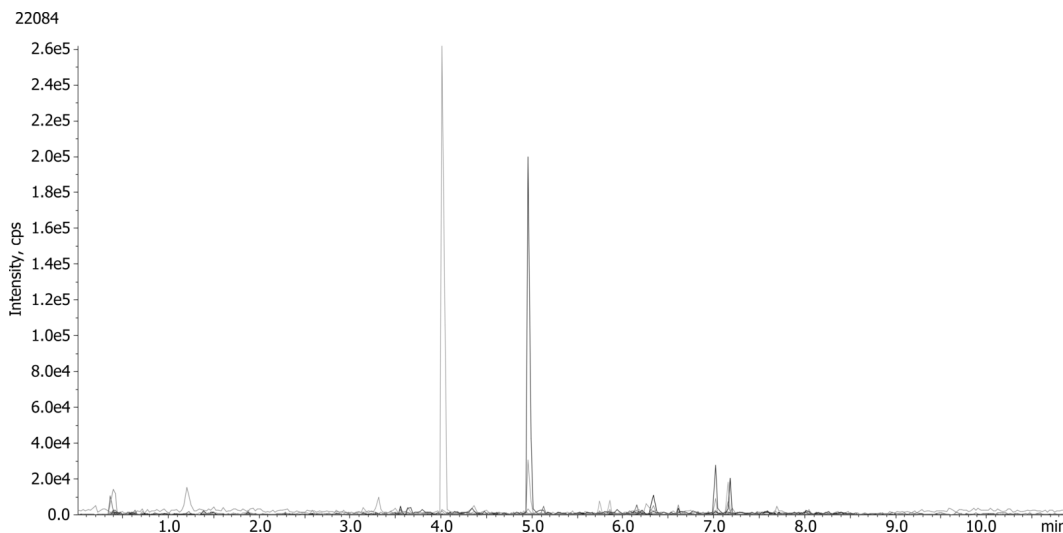


ns spiked Ground Meat Extract for antibiotic analysis on Kinetex XB-C18**Column:** Kinetex® 2.6 µm XB-C18 100 Å, LC Column 100 x 2.1 mm, Ea**Dimensions:** 100 x 2.1 mm ID**Order No:** 00D-4496-AN**Elution Type:** Gradient**Eluent A:** 0.1% Formic Acid in DI H₂O**Eluent B:** 0.1% Formic Acid in 50:50 MeOH:ACN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	1	95	5
	3	4	50	50
	4	6	5	95
	5	7.5	5	95
	6	7.51	95	5
	7	10	95	5

Flow Rate: 450 mL/min**Col. Temp.:** 50 °C**Detection:** Tandem Mass Spec (MS-MS) @ (1 °C)**Detector Info:** API 5000**ANALYTES:****1** N/A**Products used in this application:**