

Pain Panel Benzodiazapines using Kinetex 2.6u Biphenyl 50x3.0mm

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4622-Y0

Elution Type: Gradient

Eluent A: Water with 0.1% formic acid

Eluent B: Methanol with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	0.5	90	10
	3	2	75	25
	4	4.5	20	80
	5	4.51	15	85
	6	5.5	15	85
	7	7	90	10

Flow Rate: 0.6 mL/min


Col. Temp.: 40 °C

Detection: Tandem Mass Spec (MS-MS) @ (40 °C)

22109



Products used in this application:

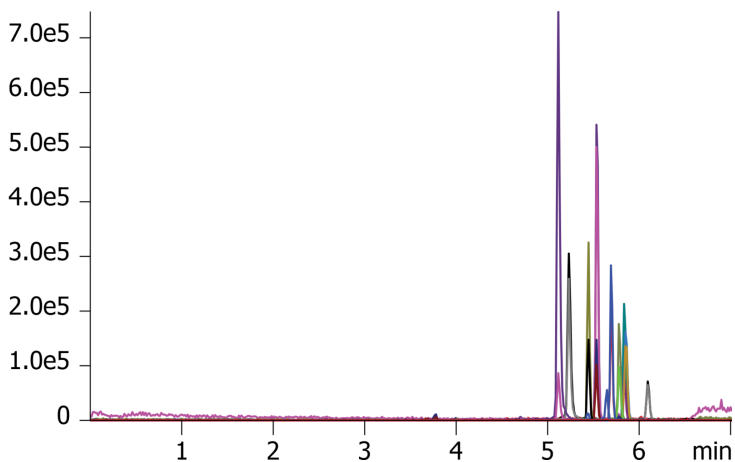

Pain Panel Benzodiazapines using Kinetex 2.6u Biphenyl 50x3.0mm

ANALYTES:

- 1** Flurazepam
- 2** Midazolam
- 3** Lorazepam
- 4** Clonazepam
- 5** Oxazepam
- 6** Hydroxyalprazolam
- 7** Nordiazepam
- 8** Flunitrazepam
- 9** Temazepam
- 10** Alprazolam
- 11** Diazepam

