

Yarra 3µm SEC-2000 Flow Rate Study PMQC STD

Column: Yarra™ 3 µm SEC-2000, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-4512-K0

Elution Type: Isocratic

Eluent A: 100 mM Sodium Phosphate pH 6.8 + 0.02% NaN₃


Gradient Profile:	Step No.	Time (min)	Pct A
1	0	100	

Flow Rate: 0.2 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 280 nm (ambient)

22125

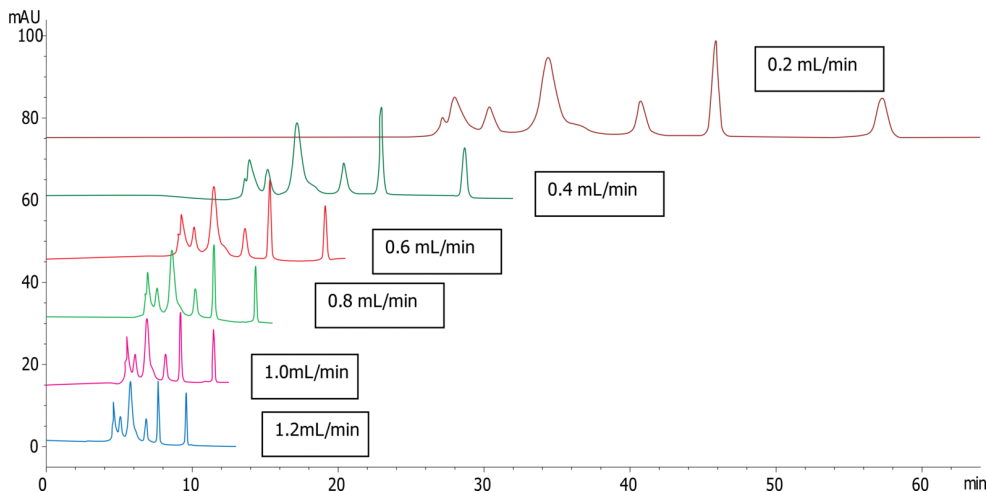

ANALYTES:

- 1 Thyroglobulin
- 2 IgA
- 3 IgG
- 4 Ovalbumin
- 5 Myoglobin
- 6 Uridine



Products used in this application:

