

Prep loading on Yarra 5u SEC-2000 PREP: Fixed concentration, increased injection volume

Column: Yarra 5 µm SEC-2000 PREP, LC Column 300 x 21.2 mm, Ea

Dimensions: 300 x 21.2 mm ID

Order No: 00H-4619-P0

Elution Type: Isocratic

Eluent A: 100 mM Sodium Phosphate pH 6.8


Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)

22128

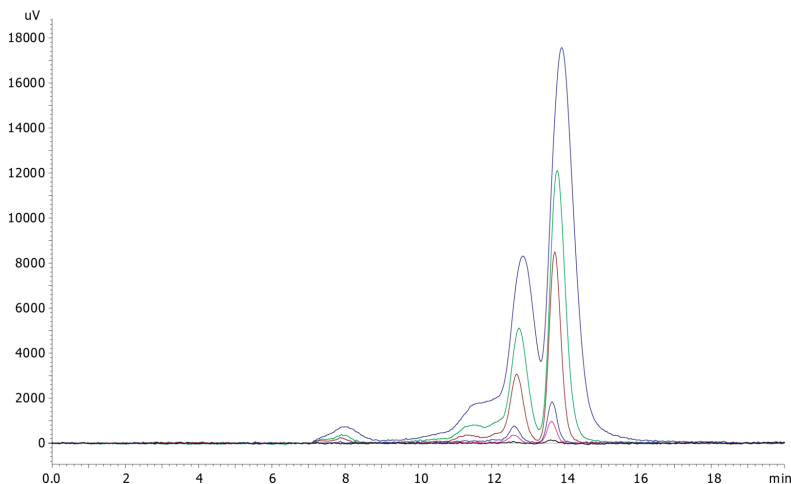

ANALYTES:

- 1 Bovine Serum Albumin
- 2 Myoglobin



Products used in this application:

