

## Aeris PEPTIDE 3.6 µm XB-C18: peptide mix, high flow rate

**Column:** Aeris™ 3.6 µm PEPTIDE XB-C18 100 Å, LC Column 50 x 2.1 mm, Ea

**Dimensions:** 50 x 2.1 mm ID

**Order No:** 00B-4507-AN

**Elution Type:** Gradient

**Eluent A:** 0.1% Formic Acid in Water

**Eluent B:** 0.1% Formic Acid in acetonitrile


| Gradient Profile: | Step No. | Time (min) | Pct A | Pct B |
|-------------------|----------|------------|-------|-------|
|                   | 1        | 0          | 95    | 5     |
|                   | 2        | 1.26       | 35    | 65    |
|                   | 3        | 1.6        | 5     | 95    |
|                   | 4        | 2          | 5     | 95    |

**Flow Rate:** 0.9 mL/min


**Col. Temp.:** ambient

**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22132

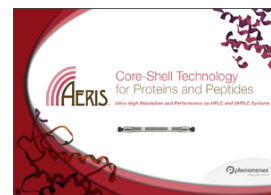

Products used in this application:

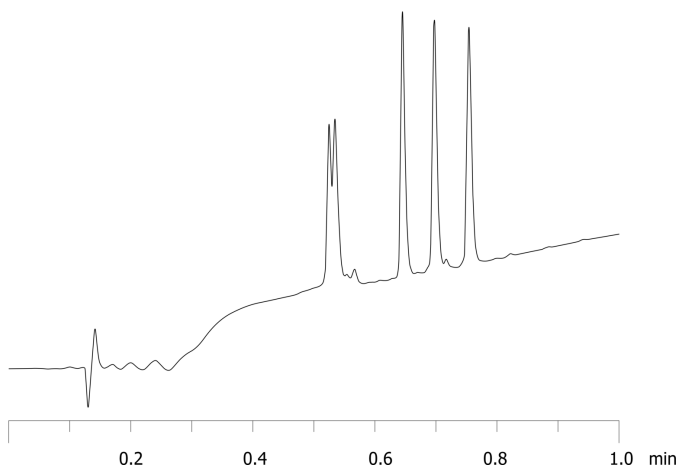

### ANALYTES:

- 1 bradykinin
- 2 dynorphin A
- 3 angiotensin II
- 4 met-enkephalin
- 5 leu-enkephalin

