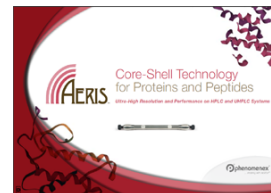
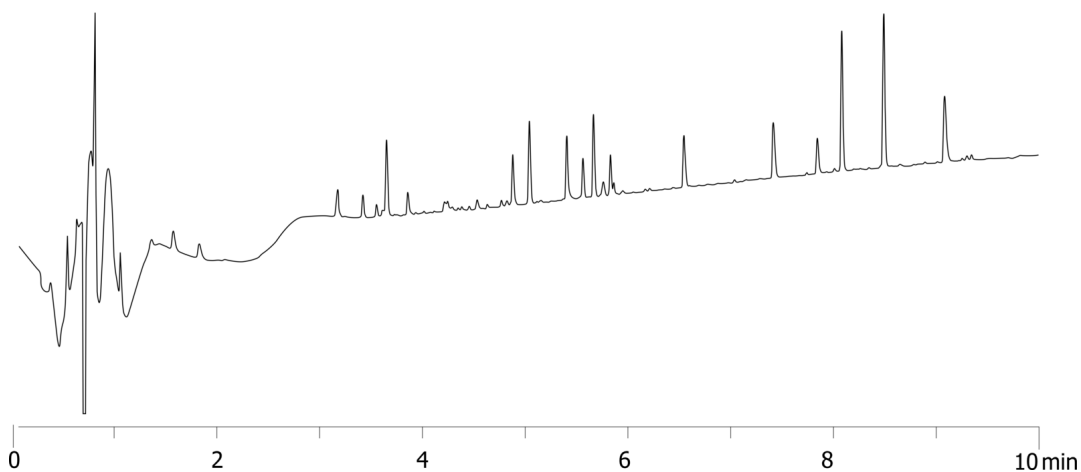


Apo-Myoglobin Tryptic Digest on Aeris PEPTIDE 1.7u XB-C18 150x2.1mm**Column:** Aeris™ 1.7 µm PEPTIDE XB-C18 100 Å, LC Column 150 x 2.1 mm, Ea**Dimensions:** 150 x 2.1 mm ID**Order No:** 00F-4506-AN**Elution Type:** Gradient**Eluent A:** 0.1% Formic Acid in Water**Eluent B:** 0.1% Formic Acid in Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	97	3
	2	1.5	97	3
	3	15	35	65
	4	17	5	95

Flow Rate: 0.5 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

22137

**Products used in this application:****ANALYTES:**

- 1 Apo Myoglobin Tryptic Digest

