

Steroid Panel on Kinetex 2.6u Biphenyl 50x3.0

Column: Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 3.0 mm, Ea

Dimensions: 50 x 3 mm ID

Order No: 00B-4622-Y0

Elution Type: Gradient

Eluent A: Water

Eluent B: Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	80	20
	2	5	40	60
	3	6	40	60


Flow Rate: 0.6 mL/min

Col. Temp.: 25 °C

Detection: UV-Vis Abs.-Variable Wave.(UV) @ (ambient)
22145



Products used in this application:

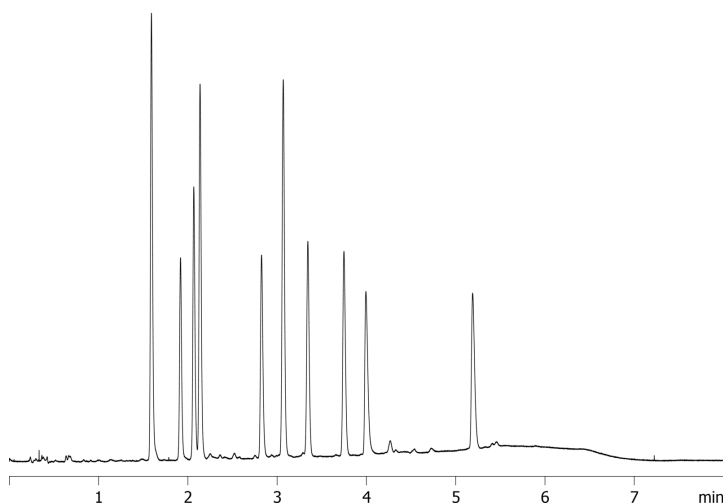

Steroid Panel on Kinetex 2.6u Biphenyl 50x3.0

ANALYTES:

- 1** Estriol
- 2** Hydrocortisone
- 3** Predisone
- 4** Cortisone
- 5** Corticosterone
- 6** Beta-Estradiol
- 7** Cortisone Acetate
- 8** 21-Hydroxyprogesterone
- 9** 17-Hydroxyprogesterone
- 10** Deoxycorticosterone

