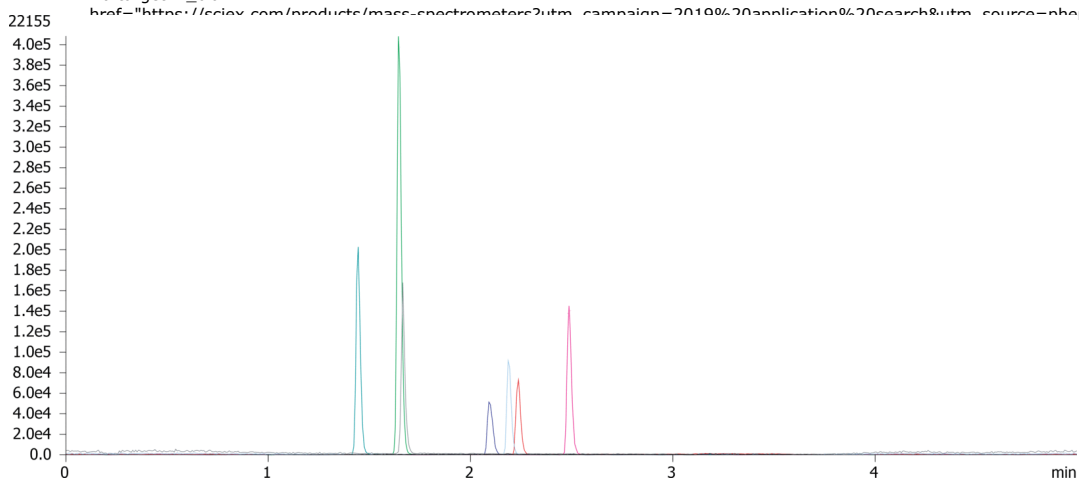


Antiepileptics in methanol using Kinetex 2.6u Biphenyl 50x2.1mm**Column:** Kinetex® 2.6 µm Biphenyl 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4622-AN**Elution Type:** Gradient**Eluent A:** Water with 0.1% formic acid**Eluent B:** Methanol with 0.1% formic acid

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	98	2
	2	3	5	95
	3	3.1	98	2
	4	6	98	2

Flow Rate: 0.5 mL/min**Col. Temp.:** 40 °C**Detection:** Tandem Mass Spec (MS-MS) @ (40 °C)**Detector Info:** SCIEX<**ANALYTES:**

- 1** Pregabalin
Retention Time: 1.44 min
- 2** Gabapentin
Retention Time: 1.65 min
- 3** Levetiracetam
Retention Time: 1.67 min
- 4** Lamotrigine
Retention Time: 2.09 min
- 5** Zonisamide
Retention Time: 2.19 min
- 6** Lacosamide
Retention Time: 2.24 min
- 7** MHD
Retention Time: 2.49 min

**Products used in this application:**