

Polystyrene MW Standard on Phenogel 5u 50A without SecurityGuard

Column: Phenogel™ 5 µm 50A, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0441-K0

Elution Type: Isocratic

Eluent A: THF

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)
22610



Products used in this application:

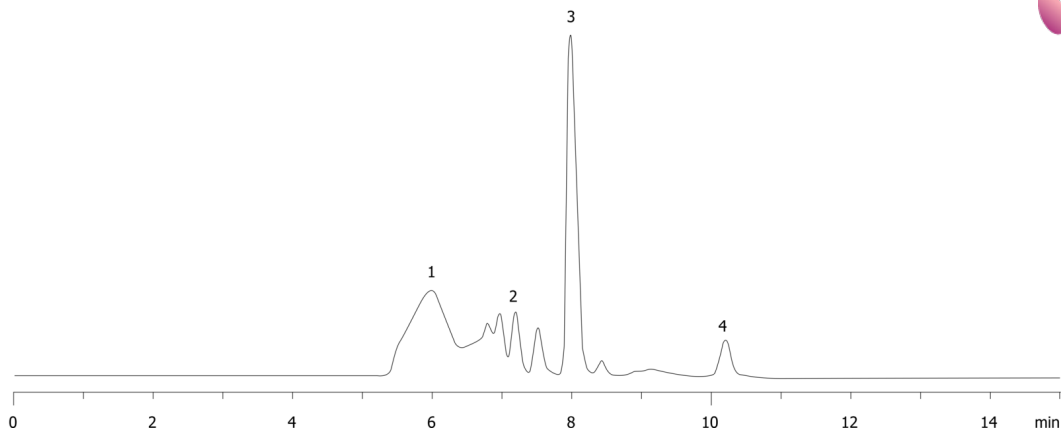

ANALYTES:

1 2400 Da

2 1490 Da

3 545 Da

4 162 Da

