


EPA 504 by GC/ECD on ZB-CLPesticides-1 (Helium Gas)**Column:** Zebtron ZB-CLPesticides-1 Column 30m x 0.32mm x 0.32 μ m , Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.32 μ m**Order No:** 7HM-G028-51**Oven Profile:** 30 °C for 2 min to 200 °C @ 30 °C/min**Carrier Gas:** Constant Flow Helium, 60 cm/sec**Injection:** Splitless :1 2 μ L @ 200°C**Detection:** Electron Capture (ECD) (220°C)**Analyst Note:** 2.0 μ L pulsed splitless @ 25 psi for 0.5 min**Products used in this application:**

22623

**ANALYTES:**

- 1 Bromodichloromethane
- 2 Chlorodibromomethane
- 3 1,2-Dibromomethane (EDB)
- 4 1,1,1,2-Tetrachloroethane
- 5 Bromoform
- 6 1,2,3-Trichloropropane
- 7 1,2-Dibromo-3-chloropropane (DBCP)

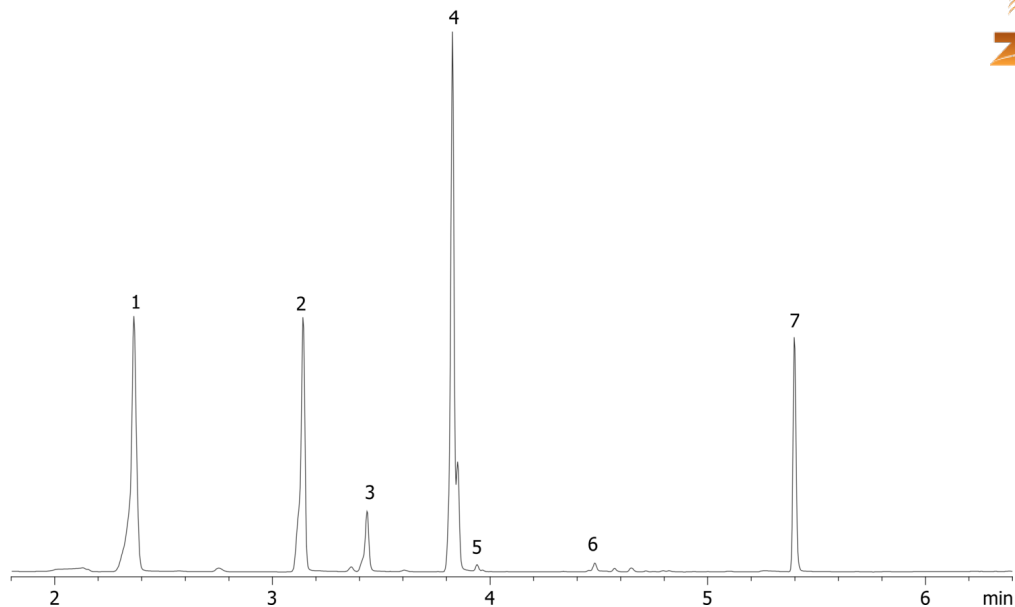