


EPA 8082 Aroclor 1242 by GC/ECD on ZB-CLPesticides-1 (Helium Gas)**Column:** Zebron ZB-CLPesticides-1 Column 30m x 0.32mm x 0.32 μ m , Ea**Phase:** Proprietary Pesticides Phase**Dimensions:** 30 meters x 0.32 mm x 0.32 μ m**Order No:** 7HM-G028-51**Oven Profile:** 120 °C to 200 @ 45 °C/min to 330 °C @ 15 °C/min for 2 min**Carrier Gas:** Constant Flow Helium, 60 cm/sec**Injection:** Splitless :1 1 μ L @ 250°C**Detection:** Electron Capture (ECD) (330°C)**Analyst Note:** 1 μ L Pulsed Splitless @ 30 psi for 20 sec.
Aroclor @ 1ppm and Surr and IS @ 100ppb in hexane**Products used in this application:****ANALYTES:**

- 1 Tetrachloro-meta-xylene [TCMX] (Surr)
- 2 Decachlorobiphenyl (IS)

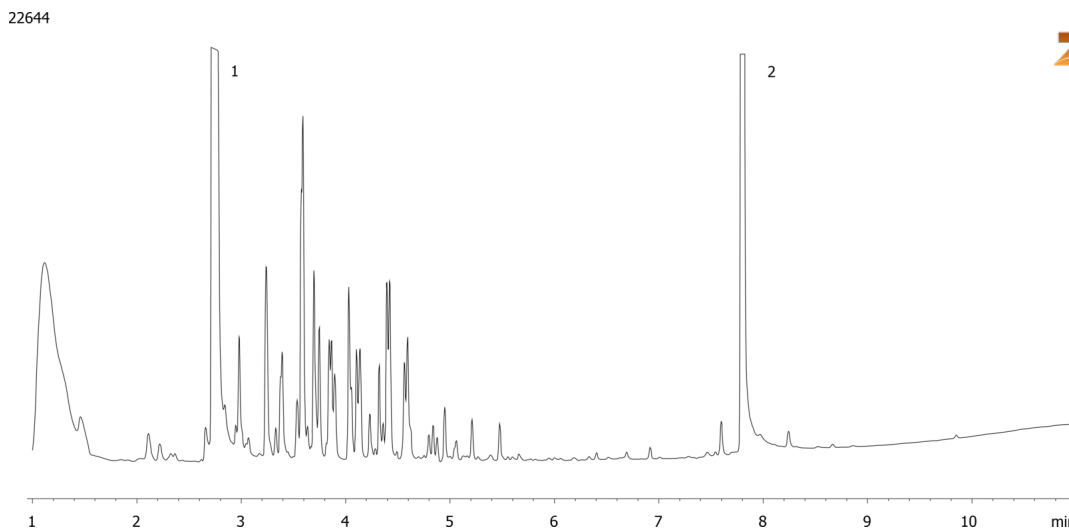