

Aeris Ig-G comparison widepore C4

Column: Aeris™ 3.6 µm WIDEPORE C4 200 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4486-E0

Elution Type: Gradient

Eluent A: 0.1% TFA


Eluent B: 0.085%TFA in AcCN

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	90	10
	2	1	80	20
	3	15	40	60
	4	15.01	90	10

Flow Rate: 1 mL/min

Col. Temp.: 40 °C


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 214 nm (ambient)

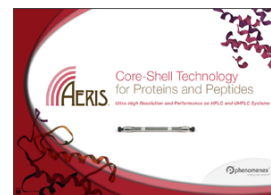

Products used in this application:



22745

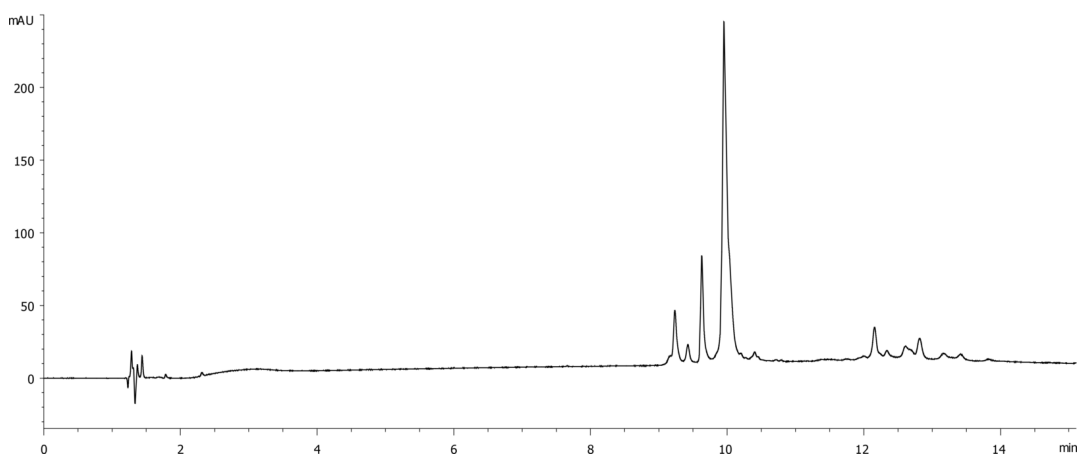

ANALYTES:

- 1 Ig-G1 Standard

