

OTC Fish Oil #1 on Phenogel 5u 500Å, 100Å, 50Å in Series

Column: Phenogel™ 5 µm 100A, LC Column 300 x 7.8 mm, Ea

Dimensions: 300 x 7.8 mm ID

Order No: 00H-0442-K0

Elution Type: Isocratic

Eluent A: THF


Gradient Profile:	Step No.	Time (min)	Pct A
	1	0	100

Flow Rate: 0.8 mL/min

Col. Temp.: 30 °C

Detection: Refractive Index (RI) @ (35 °C)

22834

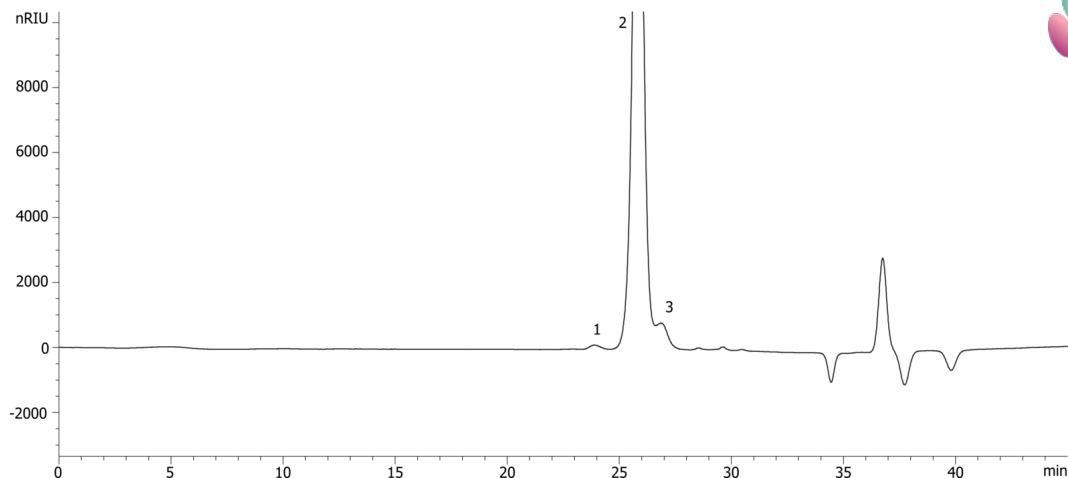

Products used in this application:



ANALYTES:

1 Oligomer

Retention Time: 23.9 min

2 OTC-1, Peak 1

Retention Time: 25.83 min

3 OTC-1, Peak 2

Retention Time: 26.87 min

