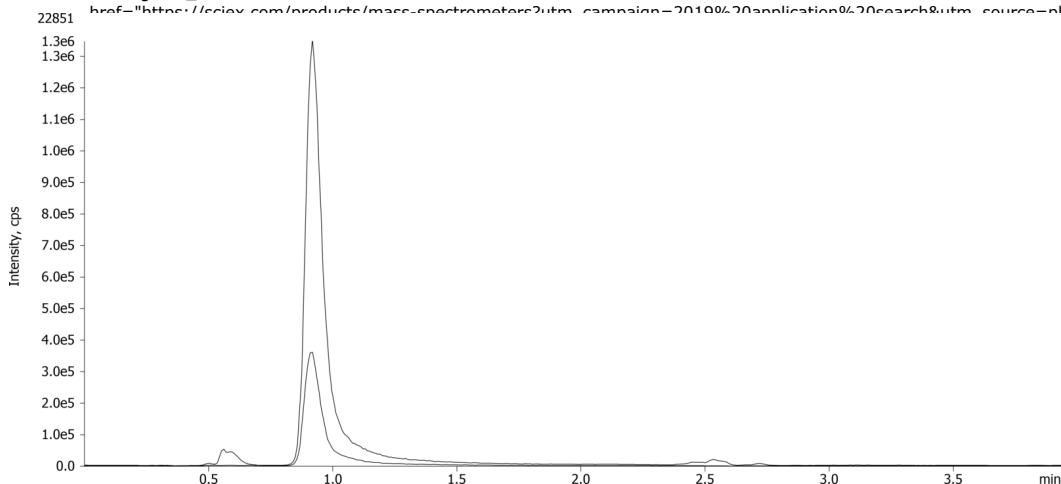


Acrylamide from Coffee by SLE and LC/MS/MS**Column:** SynergiTM 4 μ m Hydro-RP 80 Å, LC Column 50 x 2 mm, Ea**Dimensions:** 50 x 2 mm ID**Order No:** 00B-4375-B0**Elution Type:** Gradient**Eluent A:** 0.1% FA in Water**Eluent B:** 0.1%FA in MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	100	0
	2	0.6	100	0
	3	0.85	0	100
	4	3	0	100
	5	3.01	100	0
	6	5	100	0

**Flow Rate:** 0.3 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** SCIEX<**ANALYTES:****1** acrylamide

Retention Time: 0.92 min

2 acrylamide-13C3

Retention Time: 0.92 min

