

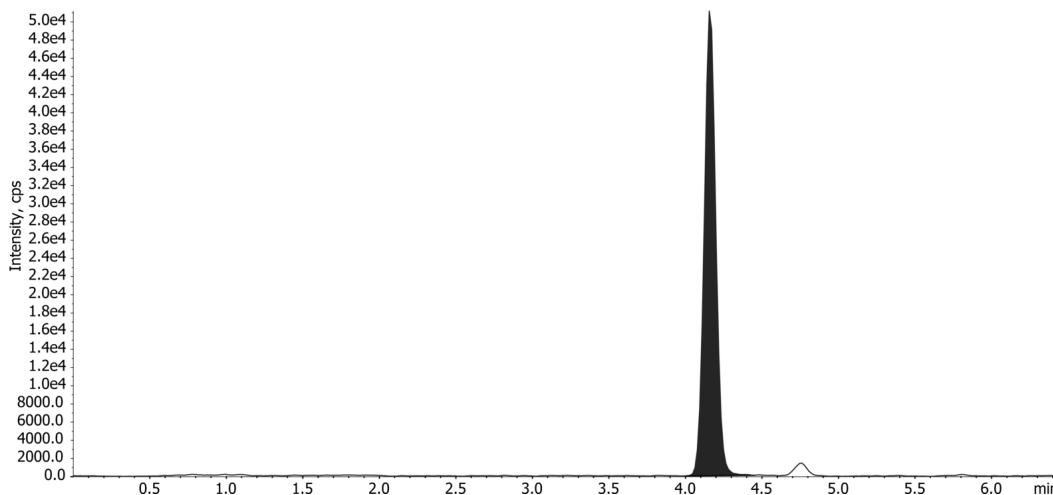
25-OH-Vitamin D2 in Human Serum using Novum SLE and Kinetex 2.6µm F5 100 x 4.6mm**Column:** Kinetex® F5 2.6µm 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4723-E0**Elution Type:** Gradient**Eluent A:** 0.1% FA in Water**Eluent B:** 0.1%FA in MeOH

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	15	85
	2	5	15	85
	3	5.01	0	100
	4	6.5	0	100
	5	6.51	15	85
	6	8	15	85

Flow Rate: 0.75 mL/min**Col. Temp.:** ambient**Detection:** Tandem Mass Spec (MS-MS) @ (ambient)**Detector Info:** <a target="_blank"

22883

href="https://sciex.com/products/mass-spectrometers?utm_campaign=2019%20application%20search&utm_source=phenomenex&utm_medium=referral">SCIEX<

**Products used in this application:****ANALYTES:****1** 25-OH-Vitamin D2

Retention Time: 4.16 min

