

Pantothenic Acid in Dietary Supplement Capsules by HPLC

Column: Synergi[™] 4 µm Hydro-RP 80 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4375-E0

Elution Type: Isocratic

Eluent A: Acetonitrile / Water / Phosphoric Acid, 85%□□3:97:0.1

Gradient


Profile:

Step No.	Time (min)	Pct A
1	10	100

Flow Rate: 1.5 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)

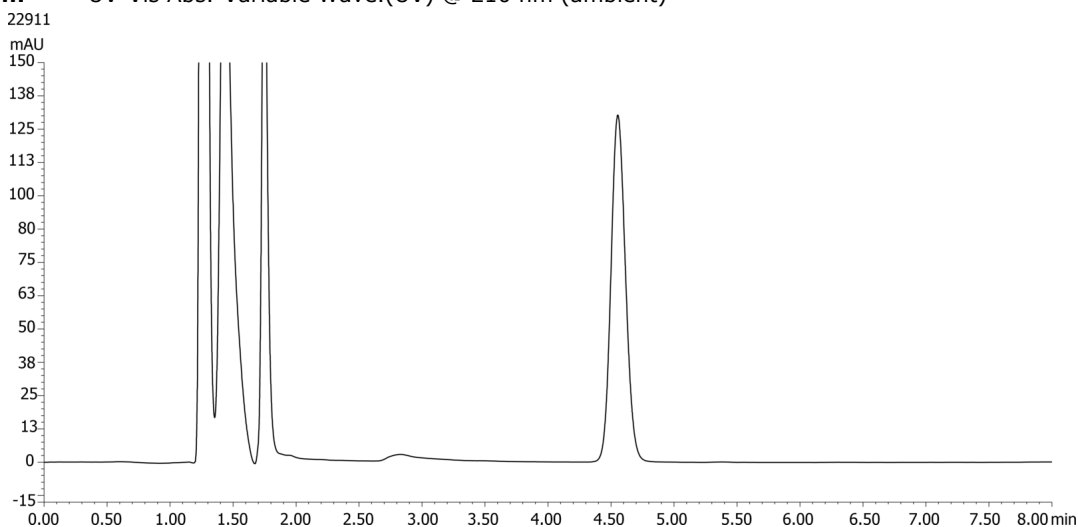

ANALYTES:

1 Pantothenic Acid

Retention Time: 4.5 min

