

Aromatic acid mix#2 on Kinetex 5u C18 150x4.6mm

Column: Kinetex® 5µm C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4601-E0

Elution Type: Isocratic

Eluent A: 0.1% H3PO4 in Water

Eluent B: Acetonitrile with 0.1% H3PO4

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	75	25


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Diode Array (PDA) @ 254 nm (ambient)



Products used in this application:

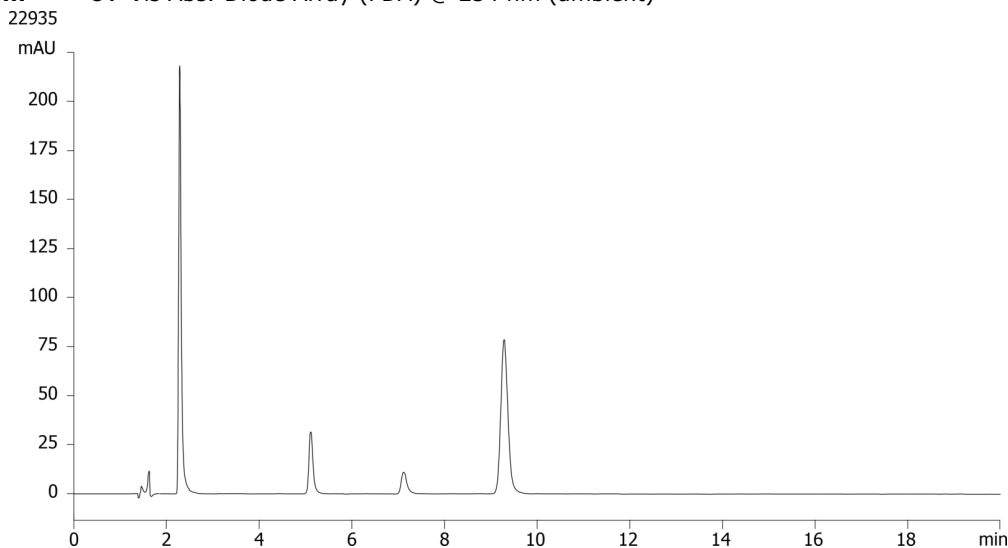

ANALYTES:

- 1 p-Hydroxybenzoic acid
- 2 Benzoic acid
- 3 o-Toluic acid
- 4 Naproxen

