

Basic drug mix on Kinetex 5u C18 150x4.6mm


Column: Kinetex® 5µm C18 100 Å, LC Column 150 x 4.6 mm, Ea
Dimensions: 150 x 4.6 mm ID
Order No: 00F-4601-E0
Elution Type: Isocratic
Eluent A: 20mM Potassium Phosphate pH 7.0
Eluent B: Methanol
Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	60	40

Profile:
Flow Rate: 0.8 mL/min
Col. Temp.: ambient
Detection: UV-Vis Abs.-Diode Array (PDA) @ 254 nm (ambient)



Products used in this application:

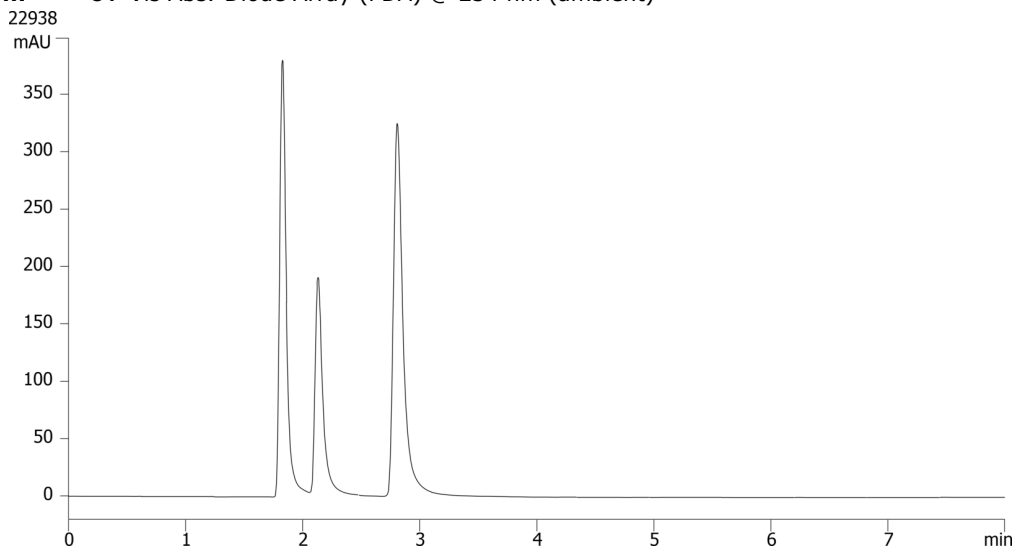

ANALYTES:

- 1 Uracil
- 2 Procainamide
- 3 N-Acetylprocainamide

