

Alpha Beta Gamma Delta Tocopherols and Tocotrienols on Kinetex 2.6u F5 150x4.6mm

Column: Kinetex® F5 2.6µm 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4723-E0

Elution Type: Gradient

Eluent A: Water

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	15	85
	2	1	15	85
	3	15	0	100
	4	17	0	100
	5	18	15	85
	6	23	15	85

Flow Rate: 1.2 mL/min

Col. Temp.: 42 °C


Detection: UV-Vis Abs.-Variable Wave.(UV) @ 290 nm (ambient)

Analyst Note: This was ran on an Agilent 1100 optimized system.

App ID 22953



Products used in this application:

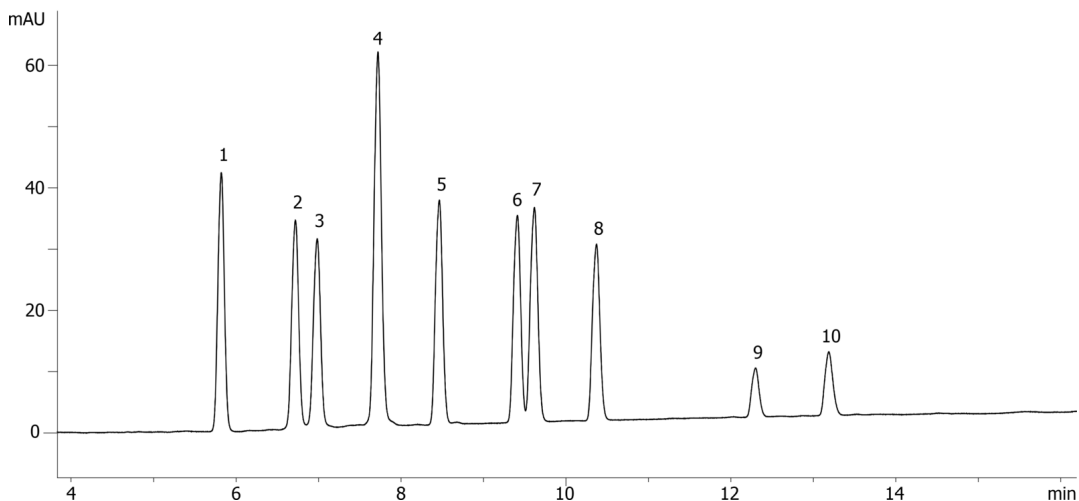

Alpha Beta Gamma Delta Tocopherols and Tocotrienols on Kinetex 2.6u F5 150x4.6mm

ANALYTES:

- 1** Delta Tocotrienol
- 2** Beta Tocotrienol
- 3** Gamma Tocotrienol
- 4** Alpha Tocotrienol
- 5** Delta Tocopherol
- 6** Beta Tocopherol
- 7** Gamma Tocopherol
- 8** Alpha Tocopherol
- 9** Vitamin E Acetate
- 10** Tocopherol Nicotinate

