


**Pheniramine and Halogenated Derivatives on Kinetex 2.6u EVO C18 100x4.6mm****Column:** Kinetex® 2.6µm EVO C18 100 Å, LC Column 100 x 4.6 mm, Ea**Dimensions:** 100 x 4.6 mm ID**Order No:** 00D-4725-E0**Elution Type:** Gradient**Eluent A:** 10mM Ammonium Bicarbonate, pH=10.2**Eluent B:** Acetonitrile

| Gradient        | Step No. | Time (min) | Pct A | Pct B |
|-----------------|----------|------------|-------|-------|
| <b>Profile:</b> | <b>1</b> | 0          | 50    | 50    |
|                 | <b>2</b> | 4          | 40    | 60    |

**Flow Rate:** 1.25 mL/min**Col. Temp.:** 25 °C**Detection:** UV-Vis Abs.-Diode Array (PDA) @ 210 nm (25 °C)**Products used in this application:****ANALYTES:**

- 1** Pheniramine
- 2** Chlorpheniramine
- 3** Brompheniramine

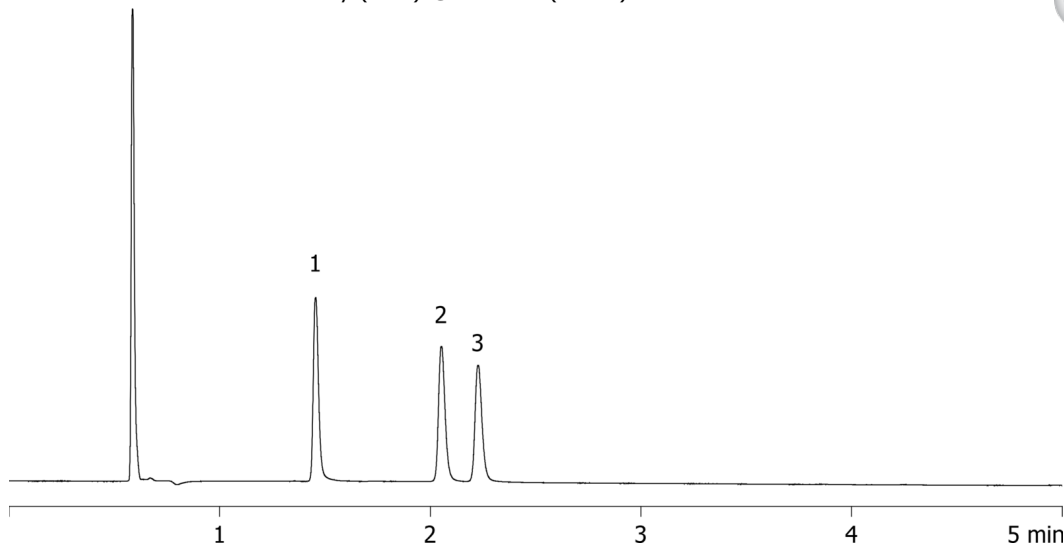