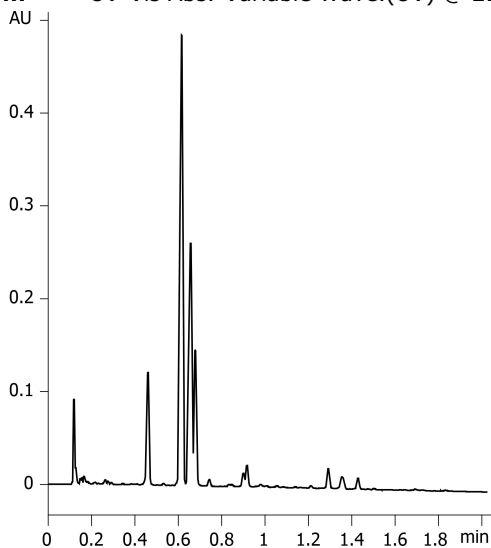


Green Coffee on Kinetex 2.6u EVO C18 50x2.1mm @ 0.8mL/min**Column:** Kinetex® 2.6µm EVO C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4725-AN**Elution Type:** Gradient**Eluent A:** 20mM Potassium Phosphate pH 2.5**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	5	30	70
	3	6	30	70
	4	6.5	95	5
	5	10	95	5

Flow Rate: 0.8 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**ANALYTES:**

- 1 Chlorogenic Acid

**Products used in this application:**