

Fast Analysis of Sugars in Fruit Juices - Orange Juice 2

Column: Rezex™ RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

Dimensions: 100 x 7.8 mm ID

Order No: 00D-0130-K0

Elution Type: Isocratic

Eluent A: Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

Flow Rate: 0.4 mL/min

Col. Temp.: 80 °C


Detection: Refractive Index (RI) @ (35 °C)
23058

REZEX™
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

REZEX™

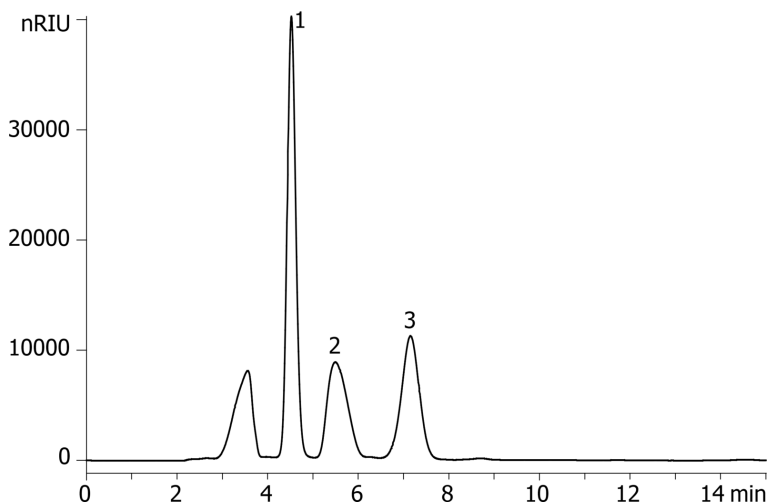

ANALYTES:

1 Orange Juice 2

