

## Fast Analysis of Sugars in Fruit Juices - Fruit Punch

**Column:** Rezex™ RCM-Monosaccharide Ca+2 (8%), LC Column 100 x 7.8 mm, Ea

**Dimensions:** 100 x 7.8 mm ID

**Order No:** 00D-0130-K0

**Elution Type:** Isocratic

**Eluent A:** Water

Gradient	Step No.	Time (min)	Pct A
Profile:	1	0	100

**Flow Rate:** 0.4 mL/min

**Col. Temp.:** 80 °C


**Detection:** Refractive Index (RI) @ (35 °C)  
23059

**REZEX™**  
HPLC COLUMNS for Carbohydrates, Oligosaccharides, and Organic Acid Analysis



Products used in this application:

**REZEX™**

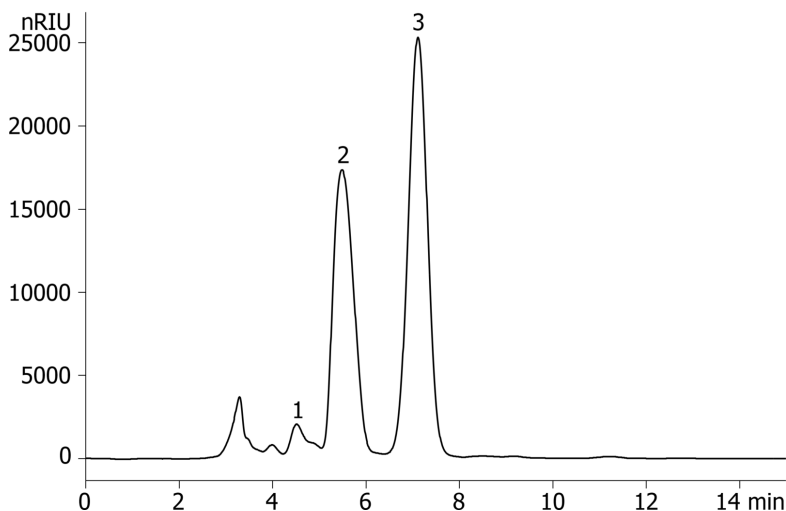

### ANALYTES:

1 Fruit Punch

