


Fluvoxamine analysis at high pH using Kinetex 1.7u EVO C18 50x2.1mm**Column:** Kinetex® 1.7µm EVO C18 100 Å, LC Column 50 x 2.1 mm, Ea**Dimensions:** 50 x 2.1 mm ID**Order No:** 00B-4726-AN**Elution Type:** Gradient**Eluent A:** 20mM Potassium phosphate pH 10**Eluent B:** Methanol

Gradient	Step No.	Time (min)	Pct A	Pct B
Profile:	1	0	95	5
	2	8	10	90

Flow Rate: 0.5 mL/min**Col. Temp.:** 30 °C**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 210 nm (ambient)**Products used in this application:****ANALYTES:****1** Fluvoxamine

Retention Time: 6.763 min

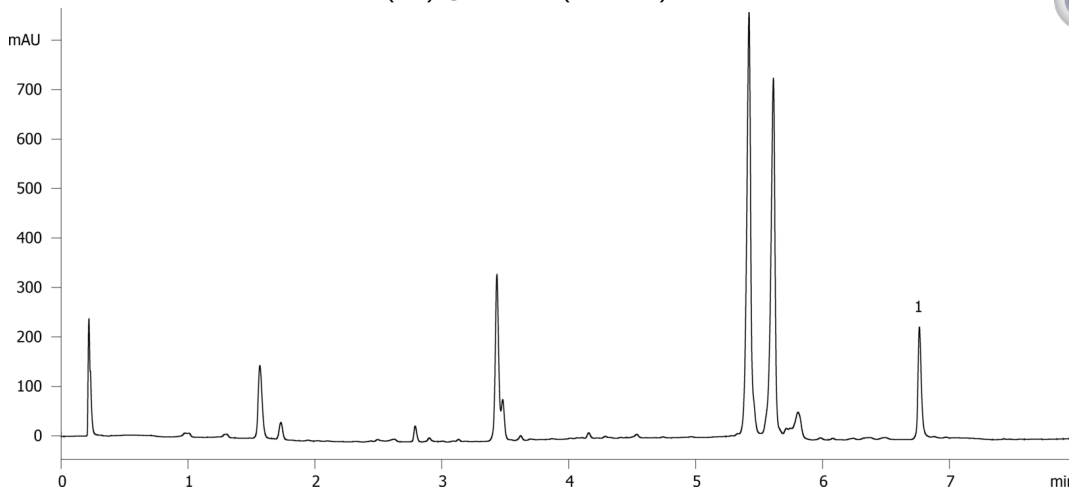