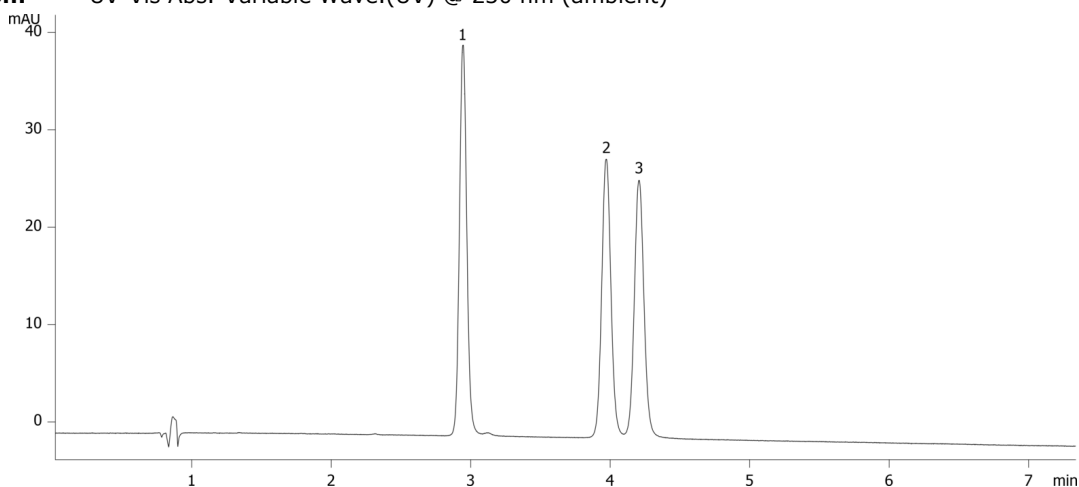


Erectile dysfunction drugs on Kinetex 5u EVO C18 150x4.6mm**Column:** Kinetex® 5µm EVO C18 100 Å, LC Column 150 x 4.6 mm, Ea**Dimensions:** 150 x 4.6 mm ID**Order No:** 00F-4633-E0**Elution Type:** Gradient**Eluent A:** 10mM Ammonium Acetate**Eluent B:** Acetonitrile

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	65	35
	2	10	50	50
	3	11	50	50
	4	12	65	35
	5	15	65	35

Flow Rate: 1.2 mL/min**Col. Temp.:** ambient**Detection:** UV-Vis Abs.-Variable Wave.(UV) @ 230 nm (ambient)**Products used in this application:****ANALYTES:****1** Tadalafil

Retention Time: 2.944 min

2 Sildenafil

Retention Time: 3.971 min

3 Vardenafil

Retention Time: 4.206 min

