

Angiotensin I,II and Aldosterone on Kinetex 5u EVO C18 150x4.6mm

Column: Kinetex® 5µm EVO C18 100 Å, LC Column 150 x 4.6 mm, Ea

Dimensions: 150 x 4.6 mm ID

Order No: 00F-4633-E0

Elution Type: Gradient

Eluent A: 5mM Ammonium Bicarbonate pH 9.0

Eluent B: Methanol

Gradient Profile:	Step No.	Time (min)	Pct A	Pct B
	1	0	95	5
	2	10	5	95


Flow Rate: 1.85 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

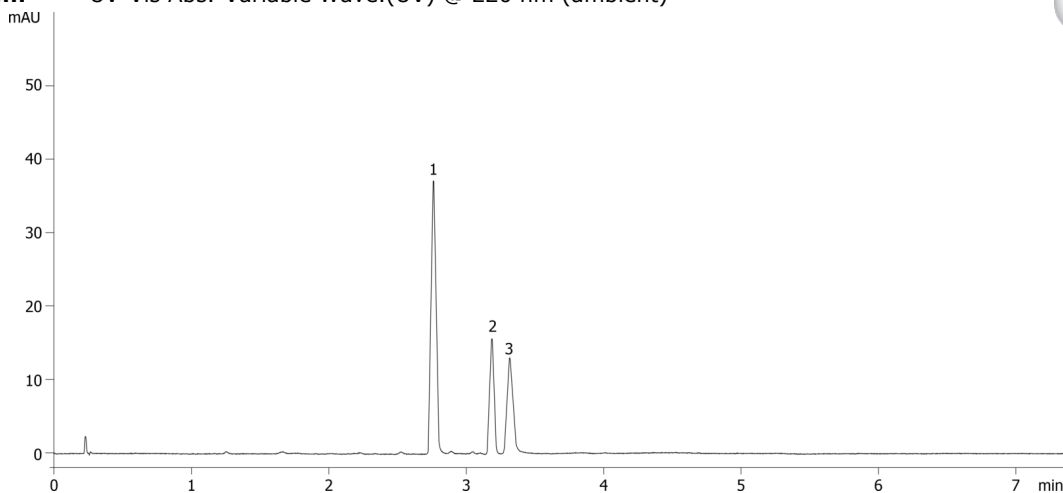

ANALYTES:

- 1 Angiotensin II
- 2 Aldosterone
- 3 Angiotensin I

