

Nadolol on LUX Amylose-1 5u 250x4.6mm in NP

Column: Lux® 5 µm Amylose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4732-E0

Elution Type: Isocratic

Eluent A: 0.1% DEA in Hexane

Eluent B: 0.1% DEA in IPA

Gradient

Step No.	Time (min)	Pct A	Pct B
1	10	60	40

Profile:


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

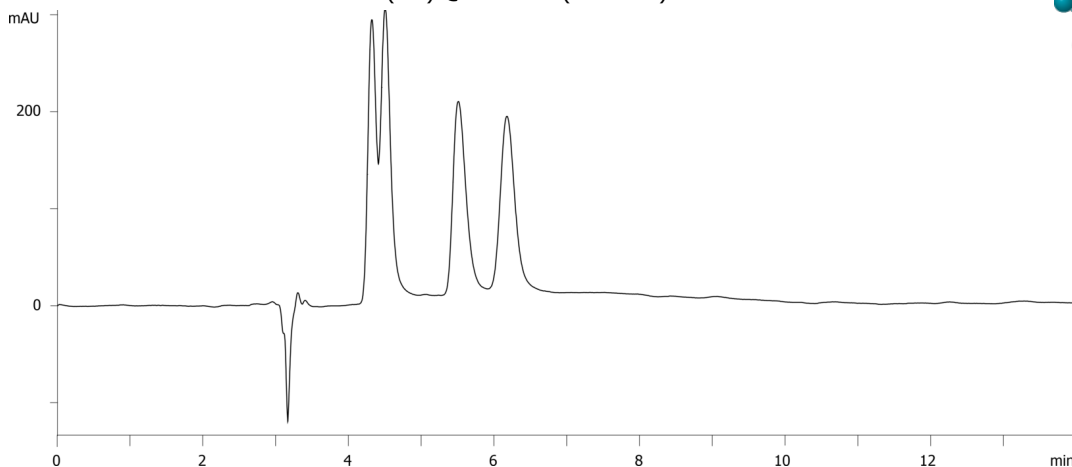

ANALYTES:

- 1 Nadolol
Retention Time: 4.313 min
- 2 Nadolol
Retention Time: 4.496 min
- 3 Nadolol
Retention Time: 5.5 min
- 4 Nadolol
Retention Time: 6.166 min

