

Pindolol on LUX Amylose-1 5u 250x4.6mm in PO

Column: Lux® 5 µm Amylose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4732-E0

Elution Type: Isocratic


Eluent A: 0.1% ethylamine in MeOH

Gradient	Step No.	Time (min)	Pct A
Profile:	1	5	100

Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 254 nm (ambient)

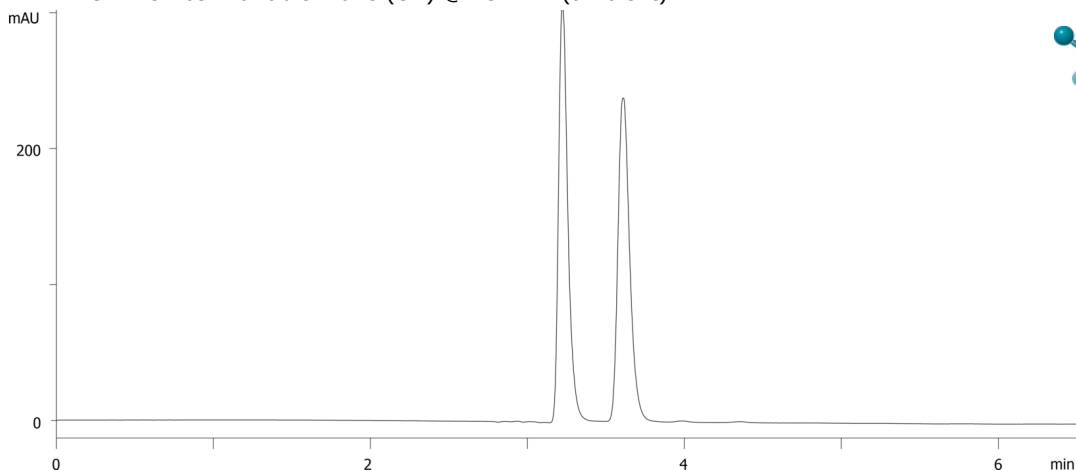

Products used in this application:



ANALYTES:

- 1 Pindolol
Retention Time: 3.225 min
- 2 Pindolol
Retention Time: 3.61 min

