

Nefopam on LUX Amylose-1 5u in RP

Column: Lux® 5 µm Amylose-1, LC Column 250 x 4.6 mm, Ea

Dimensions: 250 x 4.6 mm ID

Order No: 00G-4732-E0

Elution Type: Isocratic

Eluent A: 20mM Ammonium bicarbonate

Eluent B: Acetonitrile

Gradient

Step No.	Time (min)	Pct A	Pct B
1	0	60	40

Profile:


Flow Rate: 1 mL/min

Col. Temp.: ambient

Detection: UV-Vis Abs.-Variable Wave.(UV) @ 220 nm (ambient)



Products used in this application:

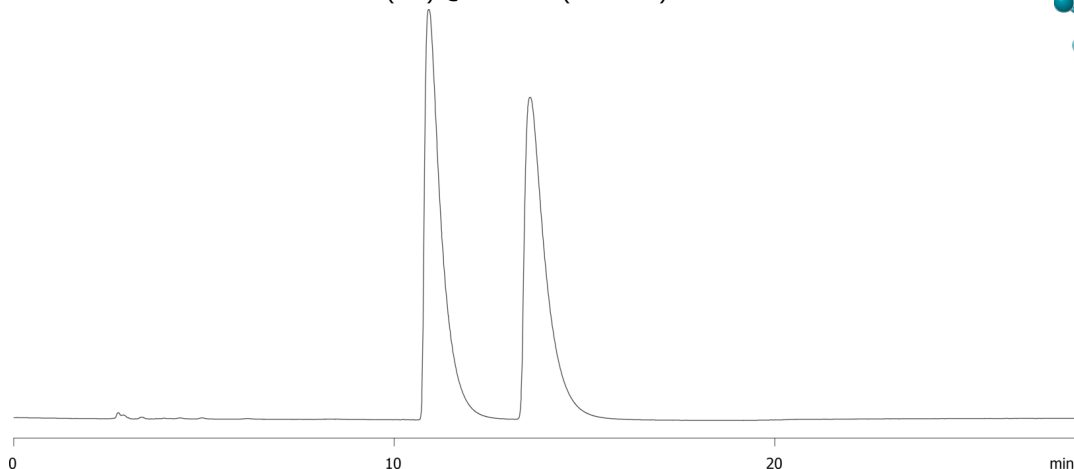

ANALYTES:

1 Nefopam

Retention Time: 10.934 min

2 Nefopam

Retention Time: 13.601 min

